

# Trauma Stretcher


Simple. Reliable. Versatile.



Stryker's Trauma Stretcher provides a secure platform for assessing and tending to critically injured and immobile patients. Its innovative design allows for **full-length X-ray capability without having to transfer or reposition** the patient.



**Slide-out cassette tray** allows you to capture precise X-ray images without moving the patient.

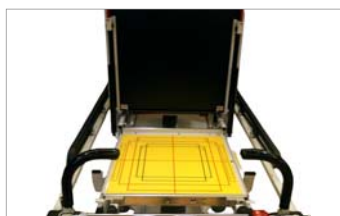

**90-degree backrest** provides the versatility to take images from a variety of positions.



Yellow



Gray



Light Blue



Red



Navy



Purple

**Multiple bumper color options** allow you to customize the stretcher for identification purposes.

## Full-length X-ray capability

Securely-mounted, lateral **C-spine cassette holder**

Precisely position **slide-out tray** for supreme imaging results

**90-degree backrest** helps facilitate imaging access

3800 E. Centre Ave.  
Portage, MI 49002 USA  
t: 269 329 2100 f: 269 329 2311  
toll free: 800 787 9537

www.stryker.com

# Trauma Stretcher



## Standard Features

- 3" (8 cm) mattress
- 8" (20 cm) Omni Surface™ casters with wheel covers
- Dual cylinder pneumatic backrest
- Dual mounted collection device hooks
- Dual pedestal hydraulics
- Dual side hydraulic controls
- Fiberesin pneumatic backrest
- Stationary foot-end litter
- Fold-down siderails
- Four IV receptacles
- Central-locking brakes with head and foot-end brake & steer pedal
- Full-length X-ray cassette tray
- Integrated storage tray with oxygen bottle holder
- Retractable fifth-wheel steering
- Rubber bumpers
- X-ray centering label
- Constant descent hydraulics

## Specifications

Model Number	1037 <b>CE</b>
Overall Length	83" (211 cm)
Overall Width (siderails up)	31.5" (80 cm)
Weight Capacity	500 lbs. (227 kg)
Height Range (to litter top)	
High	37.5" (95 cm)
Low	24.5" (62 cm)
Litter Positioning	
Backrest	0° - 90°
Trend./Reverse Trend.	±18°
Patient Surface	26" x 76" (66 x 193 cm)
Siderails	13" x 55" (33 x 140 cm)
Caster Diameter	8" (20 cm)

*Stryker reserves the right to change specifications without notice.*

## Optional Features

- Big Wheel® Advanced Mobility
- 2" (8 cm) mattress
- 4" (10 cm) Ultra Comfort Mattress
- 5" (13 cm) Pioneer® Pressure Redistribution Mattress
- C-spine cassette film holder
- Defibrillator tray/footboard/extender
- Footboard
- Footboard/charholder
- Backrest X-ray cassette holder
- IV Caddy
- IV poles
  - Permanent folding two- or three-stage
  - Tethered
  - Standard removable
- Restraint straps
- Serving/instrument tray
- Siderail pads
- Upright oxygen bottle holder
- Integrated patient transfer board
- X-ray tray centering guide
- X-ray grid and holder
- Extended warranty available