
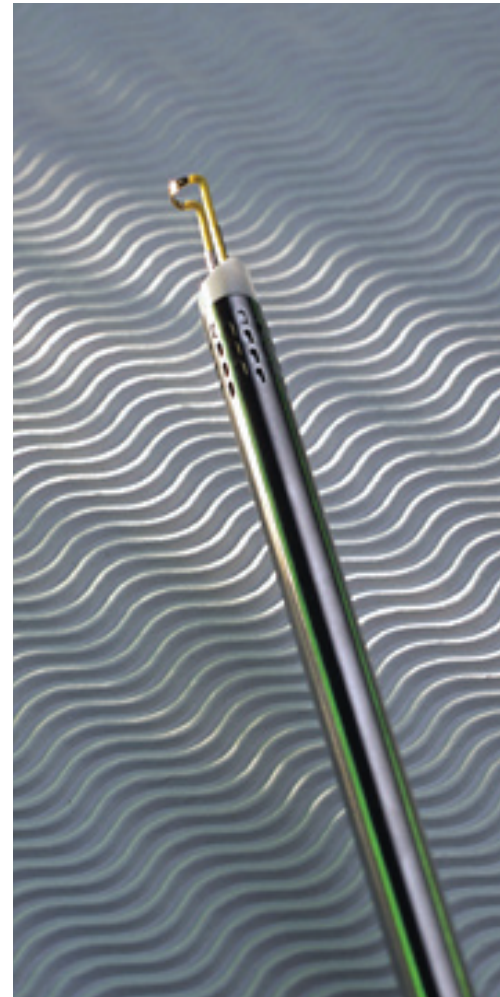



Gynecology

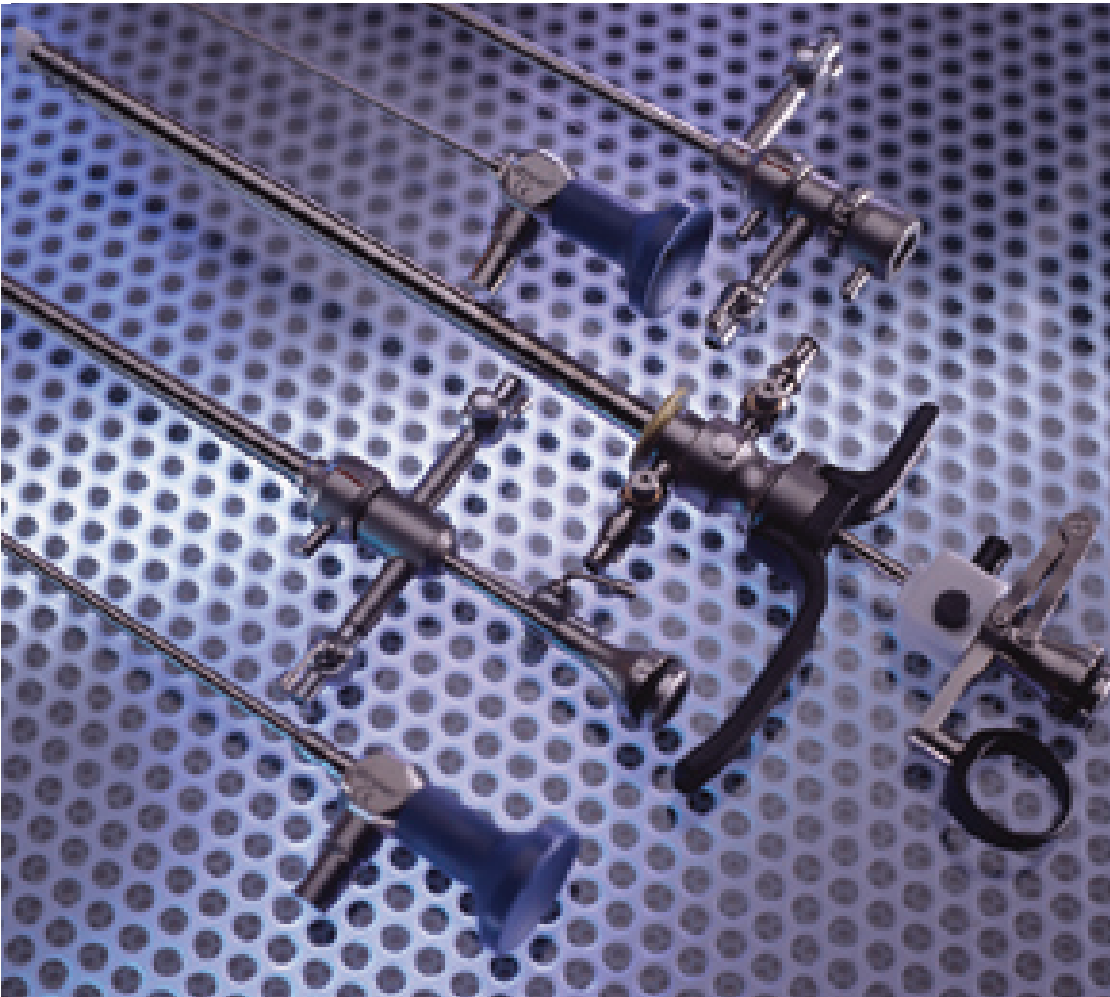

See the Difference


stryker[®]
ENDOSCOPY

Stryker's Gynecology Set is comprised of superior optics and reliable diagnostic and operative hardware; it also features the newest resectoscope on the market with high performance electrodes. Building on decades of experience, Stryker autoclavable hysteroscopes provide a bright, high-resolution image accompanied by durable and intuitive hardware. At only 188 grams, the resectoscope is the lightest on the market and has been designed with maximum flow rate and intuitive locking mechanisms. Its smooth actuation and unique electrode design allows effortless resection and rotational stability.

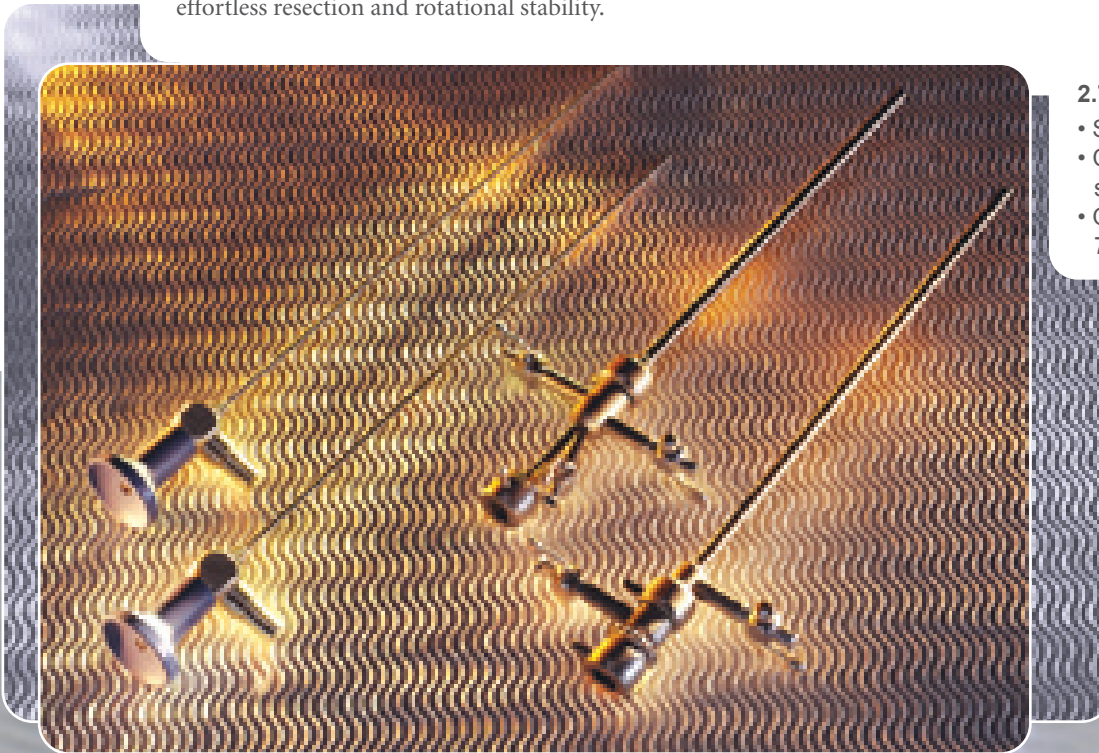


2.7mm Line

- Slide-lock connection
- Continuous flow diagnostic sheath: 5.5mm OD
- Operative sheath: 6mm OD, 7 French working channel

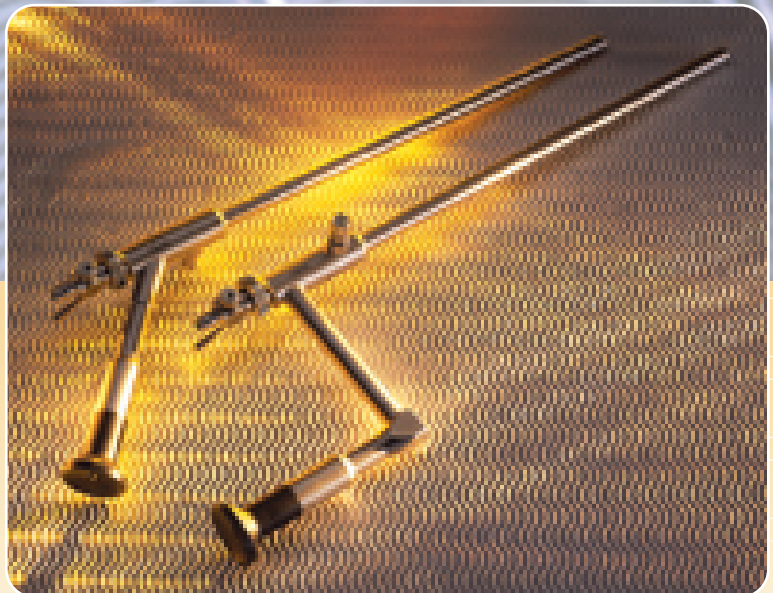


Operative Laparoscopes

- 10mm, 5mm working channel design
- 11mm, 6mm working channel design
- Durable triple tube design

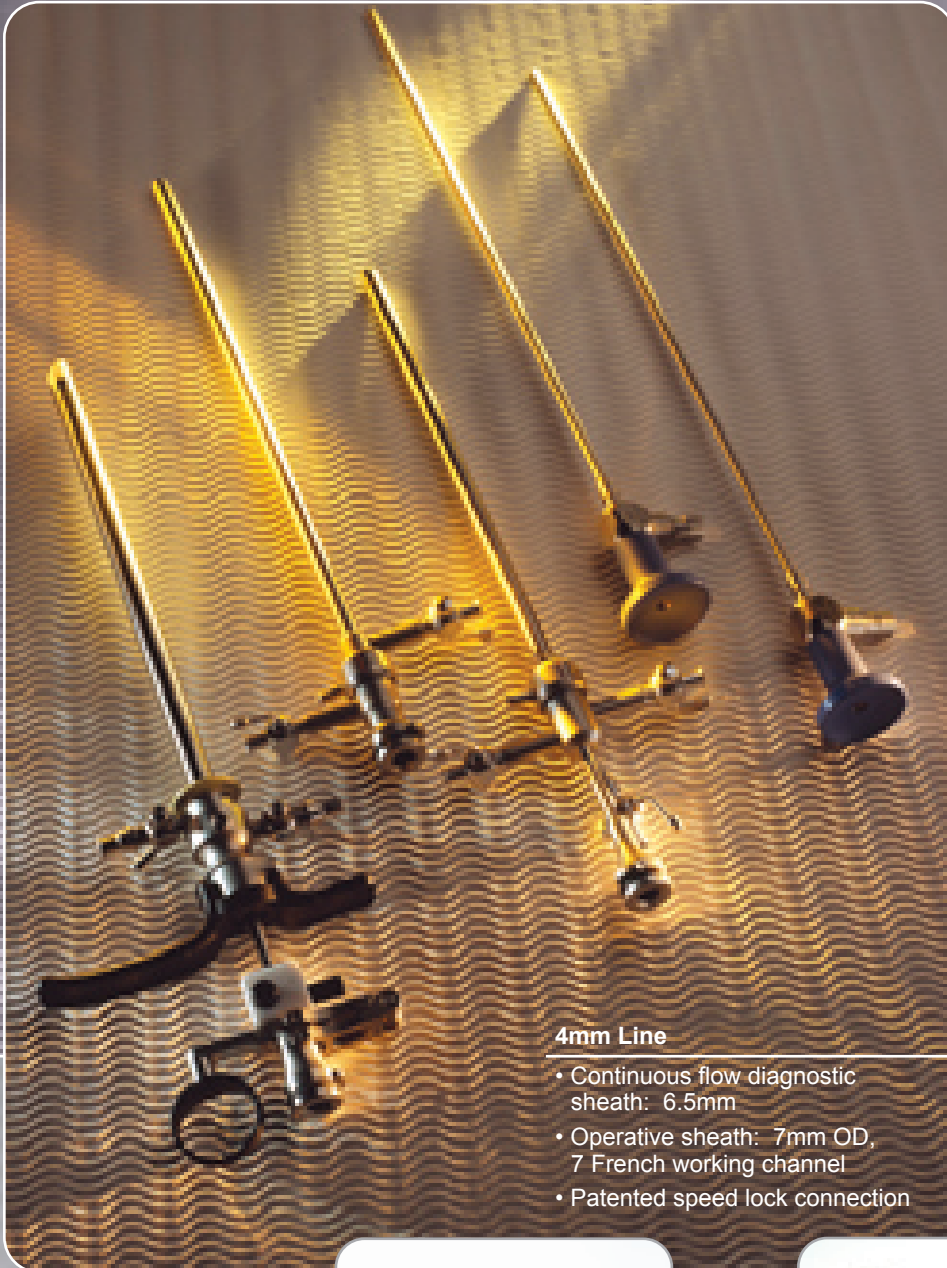

4mm Line

- Continuous flow diagnostic sheath: 6.5mm
- Operative sheath: 7mm OD, 7 French working channel
- Patented speed lock connection

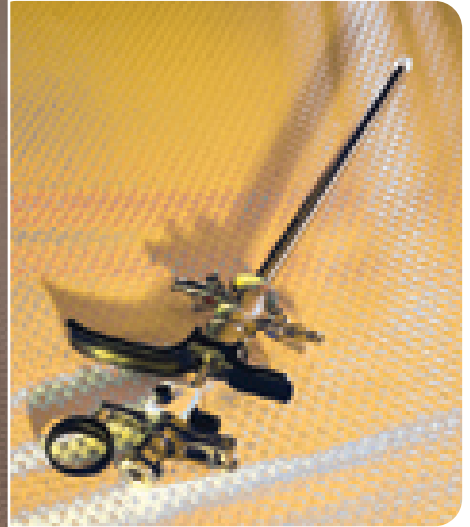


26 Fr Continuous Flow Resectoscope

- Lightweight design
- Instinctive locking mechanisms
- Smooth actuation



Instrumentation

- 7 French grasper
- 7 French punch
- 7 French scissor

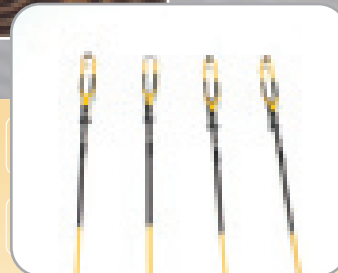

Electrodes

- Gold Resector cuts and coagulates tissue simultaneously
- Dimpled Roller Ball vaporizes tissue
- Stabilizers hug the scope, providing extra support during resection

Gynecologic Equipment and Accessories

Operative Laparoscopes

331-010-045	10mm 45° Operating Laparoscope
331-010-090	10mm 90° Operating Laparoscope
331-000-000	11mm 90° Operating Laparoscope
331-110-100	Laser Adapter
331-110-110	Merrimack Laser Coupler
331-110-120	Coherent Laser Coupler
331-110-130	Sharplan Laser Coupler
331-110-131	Sharplan 5 Start Coupler

2.7mm Diagnostic and Operative Hardware

502-743-010	2.7mm 0° Hysteroscope
502-743-030	2.7mm 30° Hysteroscope
502-944-030	2.7mm Continuous Flow Diagnostic Sheath
502-943-030	2.7mm Operative Sheath
502-958-030	2.7mm Obturator

4mm Diagnostic and Operative Hardware

502-777-012	4mm 12° Hysteroscope
502-777-030	4mm 30° Hysteroscope
502-977-030	4mm Continuous Flow Diagnostic Sheath
502-978-030	4mm Operative Sheath
502-877-030	4mm Obturator

4mm Resection Hardware

502-880-401	Working Element for Resectoscope
502-880-426	Inner/ Outer Sheath for Resectoscope
502-880-001	Standard Obturator for Resectoscope

Resectoscope Cables and Accessories

502-990-100	Aesculap/ Berchtold ESU Cable
502-990-200	Storz/ Erbe ESU Cable
502-990-300	Universal Valley Lab/ Bovis ESU Cable
233-032-881	Resection Two-Tier Tray

Resectoscope Electrodes (Sterile, box of 6)

Fit 502-880-401 working elements only

504-880-100	24 Fr Gold Coagulating Cutting Loop Electrode
504-880-200	24 Fr Roller Ball Electrode
504-880-400	24 Fr 30° Cutting Loop Electrode
504-880-600	24 Fr Dimpled Vaporizing Electrode

7 French Instruments

502-770-100	Grasper
502-770-200	Punch
502-770-300	Scissors