### *s*tryker

# Mistral-Air®

# Forced Air Warming System

### Simple. Safe. Efficient.

The Mistral-Air Forced Air System enables simple, safe and efficient management of patient temperature. The high volume blower featuring HEPA filtration is designed to work in tandem with low-pressure blankets, evenly distributing diffused warm air over the surface of the patient to help control against the onset of hypothermia.


Mistral-Air Warming System Blankets are available in numerous configurations for patient- and position-specific warming. A broad range of standard and premium reflective warming blankets facilitates use in a wide variety of clinical settings.





#### **HEPA Filter**

The Mistral-Air HEPA filter is proven as 99.99% effective in capturing what is considered to be the most difficult particle size to catch, .3 microns.

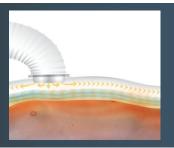

**Diffusion Technology** eliminates individual high-pressure jets of air that can cause the blanket to loft.

The blanket stays in position, keeping warm air on the patient, minimizing the disruption of air flow in the surrounding environment.



Reflective technology can reduce heat loss by up to 80% in the perioperative environment.<sup>1</sup>

Using reflective technology can help decrease patient's heat loss during transfer periods when the blanket is not connected to a blower.

#### **System Features**

#### **HEPA Filter**

Mistral-Air possesses a HEPA filter to ensure effective filtration in protective environments.

#### **Diffusion Technology**

Diffusion Technology minimizes the disruption of airflow in the surrounding environment.

#### **Reflective Technology**

Reflective Technology allows for passive warming to optimize the patient warming process during transfer.

#### Additional Features

- Filter maintenance alert
- Default startup temperature of 38° C
- Over and under temperature alarms
- Quiet, high-volume blower
- Wide range of standard and premium disposables
- Sterile blanket options





## Forced Air Warming System

*s*tryker

3800 E. Centre Ave. Portage, MI 49002 U.S.A. **t:** 269 389 2100 **f:** 269 329 2311 toll free: 800 787 9537

www.patienthandling.stryker.com

#### **Unit Specifications**

| •                              |                                   |
|--------------------------------|-----------------------------------|
| Model No.                      | MA1100-PM                         |
| Length                         | 10.9" (27.6 cm)                   |
| Width                          | 15.2" (38.5 cm)                   |
| Height                         | 9.4" (23.9 cm)                    |
| Settings                       | Low 32° C; Medium 38°; High 43° C |
| Volts                          | 110-120 V 60 Hz 6 amps            |
| Safety                         | UL60601-1, CAN/CSA C22.2 No 601.1 |
| EMC                            | EN60601-1-2                       |
| Warranty                       | 2 years                           |
| Hose                           | 6' (1.8 m)                        |
| Power Cord                     | 13' (4 m)                         |
| Mistral-Air Stand (Incl. Base) | THC5                              |

### **Warming Blankets and Suits Specifications**

#### Warming Blankets - Intraoperative

| -         | •                               |                                   |
|-----------|---------------------------------|-----------------------------------|
| Item No.  | Description                     | Size                              |
| MA2250-PM | Lower Body Plus                 | 50" x 40.2" (127 cm x 102.1 cm)   |
| MA3350-PM | Lower Body Premium Reflective   | 55.9" x 40.2" (142 cm x 102.1 cm) |
| MA2260-PM | Upper Body Plus                 | 79.1" x 29.9" (201 cm x 75.9 cm)  |
| MA3360-PM | Upper Body Premium Reflective   | 78.7" x 29.9" (200 cm x 75.9 cm)  |
| MA2265-PM | Half Upper Body Plus            | 63" x 29.9" (160 cm x 75.9 cm)    |
| MA0365-PM | Half Upper Body Reflective      | 63" x 29.9" (160 cm x 75.9 cm)    |
| MA2270-PM | Torso Plus                      | 29.9" x 48.8" (75.9 cm x 124 cm)  |
| MA2280-PM | Surgical Access (sterile)       | 86.6" x 50" (220 cm x 127 cm)     |
| MA2285-PM | Lower Surgical Access (sterile) | 65" x 50" (165.1 cm x 127 cm)     |
| MA2286-PM | Cardiac Blanket                 | 55.8" x 50" (141.7 cm x 127 cm)   |
| MA0290-PM | Surgical Access                 | 86.6" x 50" (220 cm x 127 cm)     |
| MA0510-PM | Tube                            | 68.9" x 5.9" (175 cm x 15 cm)     |
|           |                                 |                                   |

#### Warming Blankets - Underbody

| -         | -                         |                                     |
|-----------|---------------------------|-------------------------------------|
| Item No.  | Description               | Size                                |
| MA2475-PM | Pediatric Underbody Plus  | 65.9" x 35.6" (167.5 cm x 90.5 cm)  |
| MA2450-PM | Half Underbody Plus       | 66.1" x 35.4" (168 cm x 90 cm)      |
| MA3450-PM | Half Underbody Reflective | 64.2" x 40.2" (163.1 cm x 102.1 cm) |
| MA2400-PM | Full Underbody Plus       | 86.6" x 35.4" (220 cm x 90 cm)      |
| MA3400-PM | Full Underbody Reflective | 84.7" x 40.2" (215.1 cm x 102.1 cm) |
| MA2240-PM | Neonatal Plus             | 50" x 25.2" (127 cm x 64 cm)        |
| MA0340-PM | Neonatal Underbody        | 33.5" x 29.9" (85.1 cm x 76 cm)     |

#### **Warming Blankets - PACU**

| Item No.  | Description                  | Size                              |
|-----------|------------------------------|-----------------------------------|
| MA2220-PM | Adult Plus                   | 85.8" x 50" (217.9 cm x 127 cm)   |
| MA3320-PM | Adult Premium Reflective     | 85.8" x 50" (217.9 cm x 127 cm)   |
| MA2230-PM | Pediatric                    | 57.1" x 40.2" (145 cm x 102.1 cm) |
| MA3330-PM | Pediatric Premium Reflective | 57.1" x 40.2" (145 cm x 102.1 cm) |

### **Warming Suits**

| •         |             |                              |
|-----------|-------------|------------------------------|
| Item No.  | Description | Size                         |
| MA1710-PM | Small       | ≤ 65" (165 cm)               |
| MA1720-PM | Medium      | 65" – 70.9" (165 – 180 cm)   |
| MA1730-PM | Large       | 70.9" – 76.8" (180 – 195 cm) |
| MA1740-PM | Extra Large | ≥ 76.8" (195 cm)             |
|           |             |                              |

References
1. Brauer, M.D., Ph.D, D.E.A.A., A., Peri, M.D., T., English, M.D., FR.C.A., M., & Quintel, M.D., Ph.D., M. (n.d.). Perioperative Thermal Insulation. Surgical Technology International, XVI, 41-45.

\*INDICATIONS FOR USE: The Mistral-Air Warming System is a forced air warming device comprised of a warming unit and a variety of blankets. It is intended to raise and maintain patient temperature by means of surface warming.

Mistral-Air $^{\scriptscriptstyle (\!0\!)}$  is a registered trademark product of  $\label{thm:company} \mbox{ International.}$ 

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Stryker. All other trademarks are trademarks of their respective owners or holder.