

stryker

# Stair-PRO®

stair chair



# Stair-PRO®

## stair chair

As a medic you are relentless. And in your relentless pursuit for optimal patient care, you are faced with unique challenges patient handling in stairways is a challenge you may face on any given call. The strain on your back from lifting and lowering patients down stairs can add up over time. With 1 in 4 EMS workers suffering a career-ending back injury within the first four years of employment, it is clear you need an option to reduce this risk¹.

Each patient handled in a stairway can pose a risk for caregiver injury. Stair-PRO puts the power to help reduce that risk back in your hands. Stryker's Stair-PRO stair chair – **it's a powerful thing.**



### Compact folded depth

Easy-to-fold functionality allows the chair to fold to a compact size for convenient storage in existing vehicle compartments.

### Optional accessories:



O<sub>2</sub> Holder



Foot Support



IV Clip



Vinyl head support

### Key features:

#### Locking rear lift handles

Length and placement optimized for up stair mobility.

#### Stair-TREAD System

Allows operator to control movement down stairs without lifting which may help reduce the risk of lower back disorder. The dual treads, made with Kevlar®, stand up to tough use and assist descent.

#### Dual wheel locks

Stable patient transfer.

#### Upper Control Handle


Extends to be positioned for optimal operatory control.

#### Extendable foot end lift handles

Design to specific length and location that can help support ergonomic lifting technique.

#### Four-inch front caster wheels

Enable wheelchair-like mobility and maximize maneuverability in tight spaces.

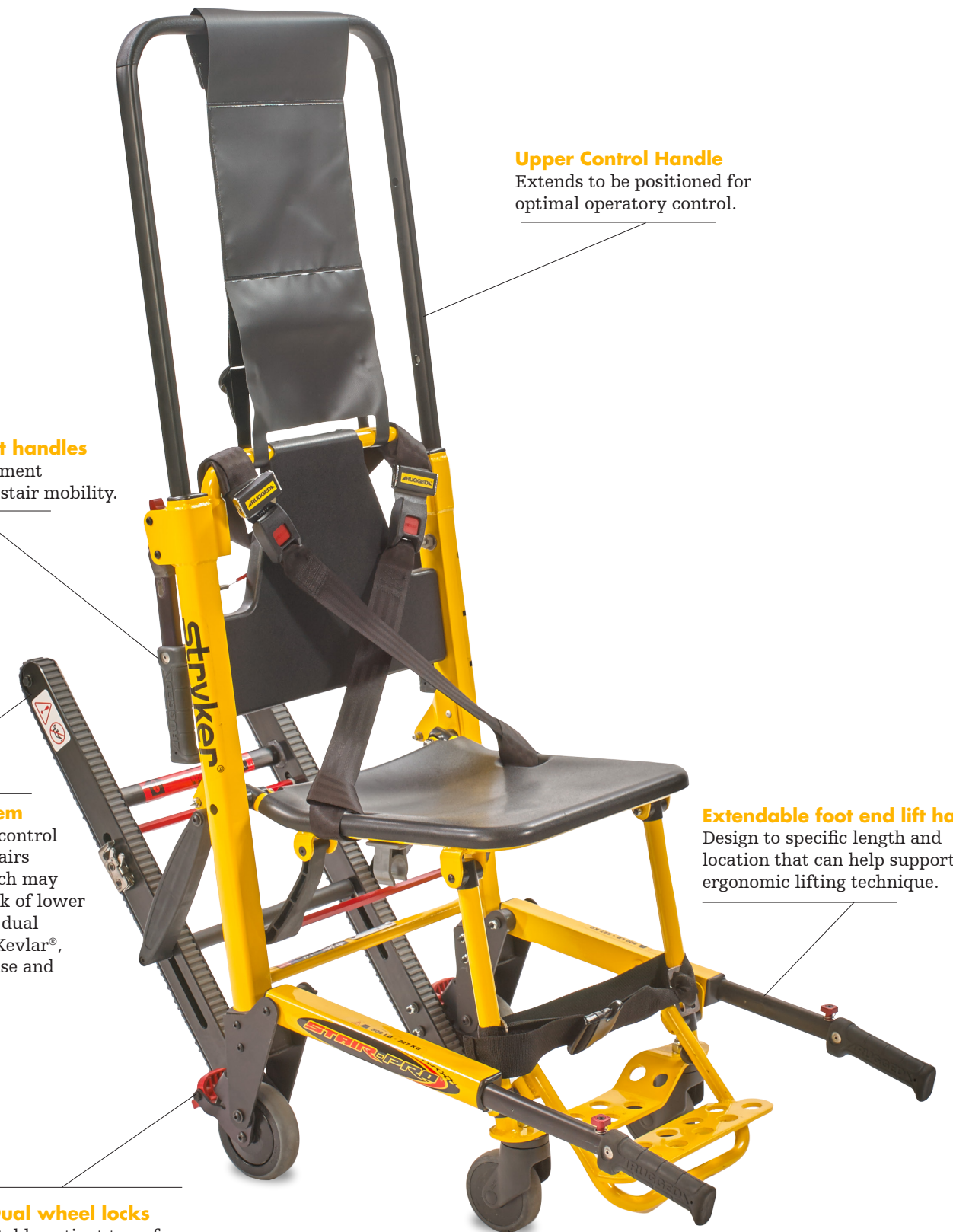

## Specifications<sup>1</sup>

• <b>Models</b>	<b>6252</b>	
• <b>Height</b>	37.5 in	(95.25 cm)
• <b>Width</b>	20.5 in	(52 cm)
• <b>Depth<sup>2</sup></b>	28 in	(71 cm)
• <b>Folded Depth</b>	8 in	(20 cm)
• <b>Weight</b>	32.8 lb	(14.9 kg)
• <b>Maximum Load<sup>3</sup></b>	500 lb	(228 kg)

<sup>1</sup> Dimensions are measured from the outermost edges of the main frame. Specifications are rounded to the nearest whole number. Conversions are calculated before rounding.

<sup>2</sup> Depth dimensions are measured with extendable handles retracted.

<sup>3</sup> Maximum load capacity is total weight distributed in accordance to basic human anatomy. EMTs must consider the weight of the patient, equipment and accessories when determining the total weight on the product.

<sup>4</sup> Patient restraint set includes two cross-chest restraints and one ankle restraint. Cross-chest restraints can also be used in a chest and lap restraint configuration.

In-service video included with every order.

Stryker reserves the right to change specifications without notice.

## References

1. Sanders, Mick J. (2011) Mosby's Paramedic Textbook (4th ed., p. 36)

## Certifications



## Warranty:

- One-year parts, labor and travel or two-year parts only
- Lifetime on all welds\*

## Extended warranties available:

- 7-year service life

## Technical Support

Stryker Technical Support comprises a team of professionals available to help with your Stair-PRO needs. Contact via phone at 1 800 STRYKER or email at [TechnicalSupport@Stryker.com](mailto:TechnicalSupport@Stryker.com).

## ProCare Service from Stryker

With maintenance, service and work-flow solutions that proactively maintain your equipment, we help you get every last penny out of your investment. Two-hour priority response callback time and 24- to 72-hour product turnaround help you focus more on your patients and less on equipment upkeep.

## Flex Financial Program

Our financial programs provide a range of smart alternatives designed to fit your organization's needs. We offer flexibility beyond a cash purchase with payment structures customizable to meet budgetary needs and help to build long-term financial stability. Specialized finance solutions for municipalities and other government agency businesses are also available.

To learn more about the wide range of financial services, call 1 888 308 3146.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: **Flex Financial, ProCare, Stair-PRO, Stryker**. All other trademarks are trademarks of their respective owners or holder.

The yellow and black color scheme is a registered trademark of Stryker Corporation

3800 E. Centre Avenue  
Portage, MI 49002 USA  
t: 269 329 2100  
toll free: 800 327 0770

[ems.stryker.com](mailto:ems.stryker.com)

Printed in U.S.A.

Mkt Lit-305 01 APR 2009 Rev D.1

Copyright © 2017 Stryker