

## ComfortGel SE

### Stretcher Support Surface



**ComfortGel SE assists in improving patient outcomes** through pressure redistribution and enhanced comfort.<sup>1</sup>


#### Helping transform the standard of stretcher support surfaces

When we developed the SE series of stretcher surfaces, we focused our design on aiding in proactive prevention.

ComfortGel SE was designed with the overall patient experience in mind and has specific features to help redistribute pressure and enhance patient comfort.

#### Key Features

- Aids in the treatment and prevention of pressure ulcers<sup>1</sup>
- Gel design affords memory and durability
- Sloped heel section assists in offloading pressure from patient's heels
- Slip-resistant bottom material eliminates the need for Velcro
- Three-sided zipper enables caregivers to more thoroughly inspect the condition of internal components

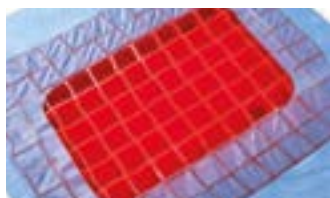

#### CoreGel Technology

Located in the sacral region, CoreGel creates a positioning pocket to help prevent patient migration.<sup>2</sup> Its design helps absorb and disperse patient weight to aid in pressure redistribution.



#### Cover design

ComfortGel SE's surface cover provides increased chemical resistance; high-frequency sealed seams, and a three-sided zipper for better visibility during internal component inspections.



#### Welded seams

Knowing surface seams are a common problem area for infection, our covers are sealed with a specific type of RF welding utilising electromagnetic bonding to limit fluid ingress and infiltration.



#### Slip-resistant bottom

ComfortGel SE's bottom features a slip-resistant material to help improve infection prevention through the elimination of Velcro while maintaining affixation to the stretcher during patient movement and transfer.



## ComfortGel SE

### Stretcher Support Surface



#### Standard Features

- CoreGel Technology
- Polycarbonate polyurethane cover
- Welded seams
- Sloped heel section
- Three-sided zipper
- Slip-resistant bottom

#### Specifications

<b>Model number</b>	1805
<b>Overall width</b>	26" (66 cm), 30" (76 cm)
<b>Overall length</b>	76" (193 cm)
<b>Thickness</b>	5.5" (14 cm)
<b>Weight capacity</b>	700 lb. (317 kg)
<b>Fire retardant inner cover*</b>	Yes

\*Flammability standards: CALTB117; 16CFR1632; CGSB CAN 2-4.2 Method 27.7-M77; BS EN 597-1; BS EN 597-2; CAL TB129; 16CFR1633; BS 6807, Clause 9; BFD IX-11; UNI9175

Stryker reserves the right to change specifications without notice.

#### Warranty

Stryker Medical, a division of Stryker Corporation ("Stryker"), warrants that its Model 1805 ComfortGel SE Support Surface Product will be free from defects in material and workmanship. Stryker's warranty covers only the following items during normal use as follows:

- Internal components: three (3) years
- Fire barrier: one (1) year
- Support surface cover: one (1) year

1. Recommended to be implemented in combination with clinical evaluation of risk factors and skin assessments made by a health care professional.
2. Best when used in combination with a stretcher equipped with Lift Assist Backrest

#### This document is intended solely for the use of healthcare professionals.

A healthcare professional must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that healthcare professionals be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A healthcare professional must always refer to the package insert, product label and/or instructions for use before using any Stryker product.

Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: ComfortGel SE, CoreGel, ShearGel, Stryker. All other trademarks are trademarks of their respective owners or holders.

The products depicted are CE marked in accordance with applicable EU Regulations and Directives.

This product is made with Intelli-Gel® hollow column configuration. "Intelli-Gel®" licensed from EdIZONE, LLC of Alpine, UT.



1362EN Rev B.0  
DFS 12/2016  
2016-12237

Copyright © 2016 Stryker  
www.stryker.com