


# LD304

## Birthing Bed

### Your partner in delivery

LD304 birthing bed provides comfort and security for both mom and caregiver.<sup>1</sup> Simple operation, mobility and an ergonomic design assist care team members while positioning to mom's comfort.

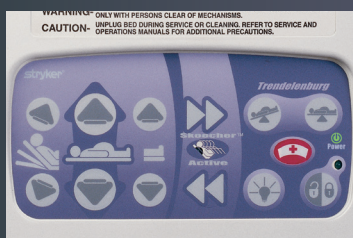

**Convenient controls**, such as single-button Trendelenburg and Stryker's powered Skoocher feature, allow customized positioning with minimal effort.



**The ergonomic design** of Stryker's foot section promotes proper body mechanics to help protect you from back strain and injury.<sup>1</sup>



**Length-adjustable footrests** provide multiple leg positioning capabilities.

**Lock-Rite foot section** glides off for quick and easy breakdown.

**Ultra Comfort Support Surface** with integrated, high-memory egg crate topper promotes comfort.

**Skoocher feature** moves mom closer to the perineal edge of the bed with the touch of a button.

# LD304

## Birthing Bed

### Standard features

- Lock-Rite slide-off foot section
- Lightweight, self-standing foot section
- Glideaway foot supports with six-way adjustability
- Single-button Trend.
- Three-position siderails
- In-rail motion controls with lockout
- 17.5-inch (44 cm) low bed height
- 5-inch (13 cm) Ultra Comfort head and seat support surface
- 3-inch (8 cm) foot mattress
- 6-inch (15 cm) Omni Surface casters with wheel covers
- Locking-caster steering
- Four-wheel, steel-ring brake system with centrally located activation pedal
- 33-inch (84 cm) mattress width with four motor functions (head, bed, foot, Trend.)
- Automatic pelvic tilt
- CPR release
- Manual backup
- Self-storing labor grips
- Plastic fluid basin
- Under-bed night light
- Wood headboard
- Permanent folding three-stage IV pole
- One year parts and labor warranty

### Specifications

<b>Model number</b>	4701
<b>Overall length</b>	93" (236 cm)
<b>Overall width</b>	
Siderails up	41" (104 cm)
Siderails down	37" (94 cm)
<b>Lift capacity</b>	
Bed	500 lb. (228 cm)
Head section	250 lb. (114 cm)
Foot section	300 lb. (136 cm)
<b>Height range (to litter top)</b>	
High	35.5" (90 cm)
Low	17.5" (44 cm)
<b>Litter positioning</b>	
Backrest	0°-70°
Pelvic tilt	4°
Maximum Trend.	8°
<b>Patient surface</b>	33" x 81" (84 x 206 cm)
<b>Support surface thickness</b>	
Head/seat	5" (13 cm)
Foot	3" (8 cm)
<b>Seat depth</b>	
Standard	12.5" (32 cm)
With Skoocher	12.5-8.5" (32-22 cm)
<b>Caster diameter</b>	
Standard	6" (15 cm)
Optional	8" (20 cm)

### Optional features

- Skoocher adjusts seat depth from 12.5 inches to 8.5 inches (32 cm to 22 cm)
- Choice of removable or permanent calf supports
- Lumbar pillow and seat pillow
- Basin liners
- Wallsaver quick-release connecting kit
- Siderail communications — includes nurse call only or nurse call with speaker, TV, radio, volume, room/read lights
- Smart TV — includes channel up/down, mute
- Pendant control
- Pendant port
- 8-inch (20 cm) casters
- Labor bar
- Linen for head and foot mattresses
- Extended warranty

### Electronics

<b>Agency approvals</b>	C-US UL (120 V AC)
<b>Volts</b>	120 VAC (Domestic) 230 VAC (International)
<b>Ampere rating</b>	10 A
<b>Current leakage</b>	< 300 microamperes
<b>Frequency</b>	60 Hz (Domestic) 50/60 Hz (International)
<b>Hospital grade plug</b>	Yes
<b>EMC Compatible</b>	Yes
<b>Motor function</b>	Bed, head, foot, Skoocher (optional), Trend.

Stryker reserves the right to change specifications without notice.

**References:** 1. Data on file. Stryker. Clinical Eval-035-B.

A healthcare professional must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. We do not dispense medical advice and recommend that healthcare professionals be trained in the use of any particular product before using the product. A healthcare professional must always refer to the package insert, product label and/or instructions for use.

Stryker or its affiliated entities own, use, or have applied for the following trademarks or service marks: Lock-Rite, Skoocher, Stryker. The absence of a product, feature, or service name, or logo from this list does not constitute a waiver of Stryker's trademark or other intellectual property rights concerning that name or logo. All other trademarks are trademarks of their respective owners or holders.

3800 E. Centre Avenue  
Portage, MI 49002 USA  
t: 269 329 2100  
toll free: 800 327 0770

**stryker.com**