

Prime X[®] X-ray Stretcher

Mobile. Safe. Efficient.

The Prime X X-ray Stretcher is the next generation of Prime Series Stretchers, loaded with all of the features you already love — now, with X-ray capability. Prime X utilizes Clearview Technology in a dual-deck design and offers a similar image quality to an x-ray table with the convenience of a mobile platform.¹



Clearview Technology

Is structural foam encased by two layers of ABS composite. This platform offers comparable image quality to x-ray tables.¹



Open cassette platform

360 degree access allows positioning of cassettes from head-end to foot-end.



Integrated scale

Provides an accurate* patient weight in both kgs. and lbs.



- 1 Advanced mobility options** are the solution to various flooring and inclines, heavy patients and long transports.
- 2 Electric lift** eliminates manual pumping, allowing for quick patient access. Simply step on the hydraulic lift pedal to raise the litter.
- 3 Easy ID** a red color scheme on the handles, siderails, bumpers, and accent rings easily identifies Prime X.
- 4 Power-washable design[™]** with 700 lb. safe working load.
- 5 Easy to raise foot section** offers enhanced cleanability through convenient access to the cassette platform

¹±3 lb (1.3 kg) for weights less than 100 lb (45 kg) and ±3% for weights greater than or equal to 100 lb (45 kg) To meet the scale accuracy claim, the patient surface must be in the flat position (Fowler backrest and gatch down), in the highest position, and the product cannot exceed 5 degrees of Trendelenburg/Reverse Trendelenburg. Scale should not be used for dosing
²1105 and 1115 models only

Prime X[®] X-ray Stretcher



Standard features

- Clearview Technology Platform
- 700 lb weight capacity
- Less than 21-inch low height (1105 and 1115 only)
- Litter-mounted oxygen bottle holder
- Power-washable design (1105 and 1115 only)
- Center-locking brake system with dual end activators
- Dual, side-mounted foot control hydraulics with uni-lower pedal
- Dual end brake/steer controls
- Glideaway Siderails
- 4 integrated IV receptacles
- +/- Trendelenburg
- 8-inch Omni Surface Casters with integrated wheel covers
- Pneumatic backrest/stationary foot
- Thermoformed ABS base
- Integrated utility tray
- Integrated transfer board

Warranty

Two years parts, labor and travel

Optional features



- Integrated scale system
- Powered lift
- Zoom stretcher with motorized drive system (1125)
- Compatible Support Surfaces:
 - o Ultra Comfort SE 30" x 3"
 - o Ultra Comfort SE 30" x 4"
- Pop-up handles
- Foot-end push handles (1105x only)
- Four-sided brake/steer controls (1105 and 1115)
- Three-sided hydraulic controls (1105 and 1115)
- Dual-ended siderail release
- IV pole options:
 - o Permanent two- or three-stage
 - o Removable
- Serving tray
- Defibrillator tray/foot extender/charting
- Chartholder
- IV caddy
- Serving tray holder/footboard
- Siderail pads
- Upright oxygen tank holder
- Restraint straps
- Retractable cord reel (1125X)

Specifications

Model number	1105X (Prime X Stretcher with Fifth Wheel)* 1115X (Prime X Stretcher with Big Wheel Advanced Mobility)* 1125X (Prime X Stretcher with Zoom Motorized Drive System)*
Overall length	
1105X and 1115X	85" (215.9 cm)
1125X	86" (218.4 cm)
Overall width	
Siderails up	38" (96.5 cm)
Siderails down (1105X and 1115X)	30.5" (77.5 cm)
Siderails down (1125X)	30.5" (77.5 cm)
Patient surface	30" x 75.25" (76 x 191 cm)
Weight capacity	700 lb. (318 kg)
Height range	
High (1105X and 1115X)	36.5" (92.7 cm)
High (1125X)	34" (86 cm)
Low (1105X and 1115X)	23.25" (59.1 cm)
Low (1125X)	23" (58 cm)
Litter positioning	
Backrest	0°–90° (±5°)
Trend./Reverse Trend. (1105X and 1115X)	±17° (±3°)
Trend./Reverse Trend. (1125X)	±16°
Siderails	14" x 57.75" (36 x 147 cm)
Caster diameter	8" (20 cm)
EMC	IEC 60601-1-2

Prime X exclusive accessories

- A. Upright x-ray cassette holder
- B. Lateral x-ray cassette holder

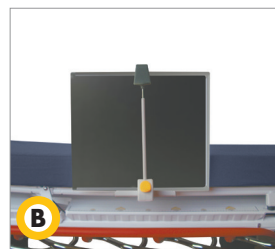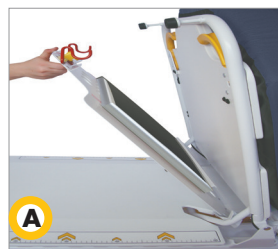

Copyright © 2017 Stryker
Mkt Lit-626 13 DEC 2011 Rev D

3800 E. Centre Avenue
Portage, MI 49002 USA
t: 269 329 2100
toll free: 800 327 0770

stryker.com

1. Prime X case study.

*Options not available in all countries. Please contact your local country representative to learn of availability in your area. Stryker reserves the right to change specifications without notice. Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: **Big Wheel, Clearview Technology, Prime X, Stryker, Zoom**. All other trademarks are trademarks of their respective owners or holder.