

# Prime TC<sup>®</sup> Transport Chair



## Patient transport redefined.


When it comes to the patient experience, details matter. That is why every aspect of the Prime TC Transport Chair is designed to help create a welcoming and safe environment that speaks volumes to caregivers, patients and their family members. Obvious touchpoints help make operation intuitive for all users. BackSmart ergonomics have been incorporated into the design to help reduce the bending and reaching that may be associated with caregiver injury. A rigid frame helps enhance durability while reducing the risk of theft. This investment in the well-being of patients and caregivers tells them you care about their comfort and safety – it tells them that they matter.



Flip-up footrests help reduce trip hazards by providing a clear path of patient ingress and egress. Foot-operated swing away functionality enhances access to the patient and closer positioning to beds, cars or tables.



BackSmart push handles help enhance caregiver safety by allowing caregivers to position their elbows at an ergonomic 90° angle and accommodate caregivers of virtually any height.

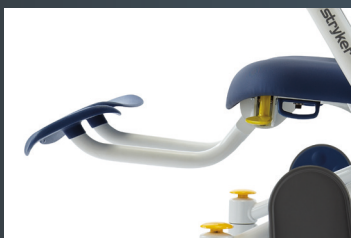

Optional leg supports are permanently attached for comfortable elevation in a fixed upright position. Easily stow away when not in use.

**Big Wheel** maneuverability helps make steering and cornering easier.

**One-touch central brake** reduces bending and reaching while enhancing caregiver efficiency.

**Stand-assist armrests** provide a raised and extended touchpoint for patient ingress and egress.

**Rigid, nestable frame** is designed to help reduce the risk of theft while preserving valuable floor space.

**Designed for cleanability**, the chair frame consists of seamless components that are easy to disinfect and are power-washable.\*

**Yellow activation points** make operation intuitive for all users.

**Movable, but not removable, components** help reduce the risk of lost or stolen parts.

**500 lb safe working load** accommodates even heavy patients.

**Durable 1" thick IV pole** is chrome-plated to help prevent rust and features enclosed hooks for safety.

**Upright O<sub>2</sub> tank holder** helps minimize bending and lifting required to place a heavy O<sub>2</sub> tank underneath the chair.

\*Chairs with the optional padded seat are not power washable.

# Prime TC Transport Chair



## Standard features

- Stand-assist armrests
- Flip-up footrests with swing-away functionality
- BackSmart push handles
- One-touch central brake
- Anti-tip wheels
- Chart holder
- Molded plastic, seamless design
- Power-washable\*
- Rigid powder-coated steel frame
- Nesting design
- Yellow activation points
- Big Wheel maneuverability
- Non-removable parts

## Optional features

- Leg supports
- Upright O<sub>2</sub> tank holder
- Chrome-plated, 1" thick IV pole
- IV pole topper available in six color options
- Foley bag hooks
- Padded seat\*
- Lap belt

\*Chairs with the optional padded seat are not power-washable

This information is intended solely for the use of healthcare professionals. A healthcare professional must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. We do not dispense medical advice and recommend that healthcare professionals be trained in the use of any particular product before using the product.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A healthcare professional must always refer to the package insert, product label and/or instructions for use. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your representative if you have questions about the availability of Stryker's products in your area. Stryker reserves the right to change specifications without notice.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Big Wheel, BackSmart, Prime TC, Stryker. All other trademarks are trademarks of their respective owners or holders. The absence of a product, feature, or service name, or logo from this list does not constitute a waiver of Stryker's trademark or other intellectual property rights concerning that name or logo.

## Specifications

|                        |                          |
|------------------------|--------------------------|
| Model number           | 1460                     |
| Overall length         | 40.2" (102.1 cm)         |
| Overall width          | 28.3" (71.8 cm)          |
| Seat width             | 21.5" (54.6 cm)          |
| Seat height            |                          |
| From footrest          | 15" (38.1 cm)            |
| From floor             | 21" (53.3 cm)            |
| Chair height           |                          |
| With IV pole           | 73" (185.4 cm)           |
| Without IV pole        | 45" (114.3 cm)           |
| Push handle grip range |                          |
| Width                  | 16"-20.5" (40.6-52.1 cm) |
| Height from floor      | 35"-45" (88.9-114.3 cm)  |
| Weight capacity        | 500 lb (226.8 kg)        |

## Warranty

Three years parts, labor and travel

3800 E. Centre Avenue  
Portage, MI 49002 USA  
t: 269 329 2100  
toll free: 800 327 0770

[stryker.com](http://stryker.com)