


Prime X Series Stretchers


Image quality. Access. Versatility.



360 degrees of access to quickly and easily place the cassette for full body imaging.

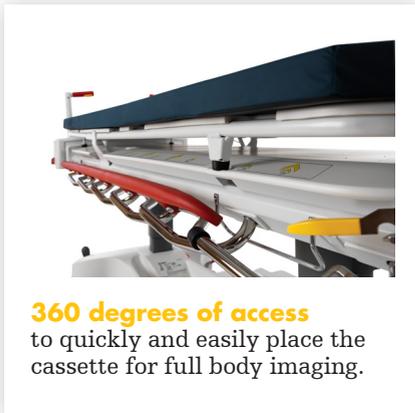

ClearView Technology combined with a dual-deck design offers the image quality you expect on a table¹ with the convenience of a mobile platform.

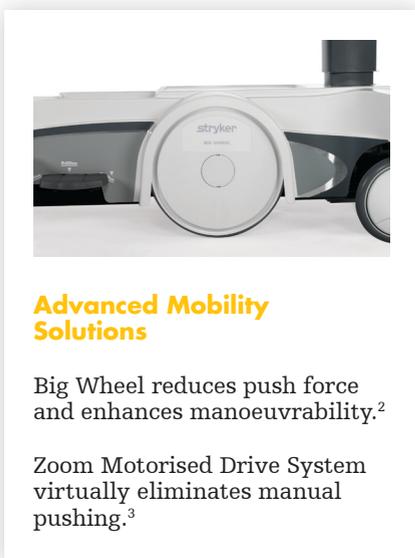

Advanced Mobility Solutions

Big Wheel reduces push force and enhances manoeuvrability.²

Zoom Motorised Drive System virtually eliminates manual pushing.³



Our Prime X offers the security and efficiency of the Prime Series stretcher, now with ClearView technology.

The dual-deck design of the patient platform allows positioning of cassettes at any point under the patient, from head to foot and side to side. Open architecture accommodates virtually any size of cassette. This versatile solution helps caregivers deliver enhanced patient care throughout the hospital by combining a mobile platform with superior image quality.

Dual-deck design allows cassette to be adjusted to the patient.

Open architecture accommodates virtually any size of cassette.

Beveled edges and visual alignment guides on the cassette platform help make cassette placement quick and easy.

Integrated scales help reduce patient transfers and make patient care efficient.

Electric lift makes raising the stretcher easy by simply stepping on the pedal.

317.5 kg (700 lb) weight capacity accommodates your heavier patients.

Articulating backrest and foot section allows complete access to the dual deck for easy cleaning.

Glideaway Siderails provide exceptional coverage and tuck neatly under the litter for unobstructed lateral imaging and storage when you don't need them.

1. Prime X Case Study. 2. Stryker Prime Series Stretcher with Big Wheel and electric lift reduces exertion when lifting the stretcher litter, 2015. 3. Zoom Motorized Drive Reduces Perceived Exertion when Transporting Up an Incline, 2013.

Prime X Series Stretchers

Image quality. Access. Versatility.

Specifications

Model number	1105X (Prime X with Fifth Wheel) 1115X (Prime X with Big Wheel) 1125X (Prime X with Zoom Drive)
Overall length	
1105X and 1115X	215.9 cm (85")
1125X	218.4 cm (86")
Overall width	
Siderails up	96.5 cm (38")
Siderails down (1105X and 1115X)	77.5 cm (30.5")
Siderails down (1125X)	78.7 cm (31")
Patient surface	76 x 191 cm (30" x 75.25")
Weight capacity	317.5 kg (700 lb.)
Height range	
High	94.6 cm (37.25")
Low (1105X and 1115X)	61 cm (24")
Low (1125X)	66.7 cm (26.25")
Litter positioning	
Backrest	0° - 90°
Foot section	0° - 40°
Trend./reverse Trend. (1105X and 1115X)	±17°
Trend./reverse Trend. (1125X)	±16°
Siderails	36 x 147 cm (14" x 57.75")
Caster diameter	20.3 cm (8")
EMC	IEC 60601-1-2

Standard features

- Thermoformed ABS base with extra capacity storage tray
- 317.5 kg. weight capacity
- 76.2 cm (30") wide patient surface with 7.6 cm (3") Enhanced Comfort mattress
- Litter-mounted oxygen bottle holder
- Integrated transfer board
- Power-washable design (1105X and 1115X only)
- Four-wheel, centre-locking brake system with dual end activators
- Dual, side-mounted foot control hydraulics with uni-lower pedal
- Glideaway siderails with integrated foot end push/pull handles
- 4 integrated IV receptacles
- +/- trendelenburg
- 20.3 cm (8") Omni Surface Casters with integrated wheel covers
- Pneumatic backrest/stationary foot
- Constant descent hydraulics
- Dual end brake/steer controls
- ClearView Technology Platform
- **Two year warranty**

Optional features

- 10.1 cm (4") Enhanced Comfort mattress
- 10.1 cm (4") or 12.7 cm (5") Ultra Comfort mattress
- Four-sided brake and steer controls (1105X & 1115X only)
- Pop-up push handles (1105X and 1115X only)
- Integrated scale system
- Three-sided hydraulic controls (1105X and 1115X only)
- Foot-end push handles (1105X only)
- Defibrillator tray
- Defibrillator tray/foot extender/charting surface
- Dual-ended siderail release
- Footboard/chartholder
- IV caddy
- IV poles (permanent, folding two- or three-stage)
- Serving tray
- Serving tray holder/footboard
- Siderail pads
- Upright oxygen bottle holder
- Restraint straps
- Retractable cord reel
- Lateral cassette
- Powered lift (1105X and 1115X only)
- Upright cassette holder
- Lateral cassette holder



1105 X • Fifth wheel



1115 X • Big wheel



1105 X • Zoom

A healthcare professional must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that healthcare professionals be trained in the use of any particular product before using it.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A healthcare professional must always refer to the package insert, product label and/or instructions for use before using any Stryker product.

Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Stryker. All other trademarks are trademarks of their respective owners or holders.

For more information please contact your local sales representative.

Stryker Australia Pty Ltd
8 Herbert Street
St Leonards NSW 2065
61 2 9467 1000
www.stryker.com.au

Stryker New Zealand Limited
PO Box 17-136
Greenlane, Auckland, New Zealand
64 9 573 1890
www.stryker.com