


stryker

SideKick[®] Disinfecting Wipes

Compatible. Versatile. Convenient.

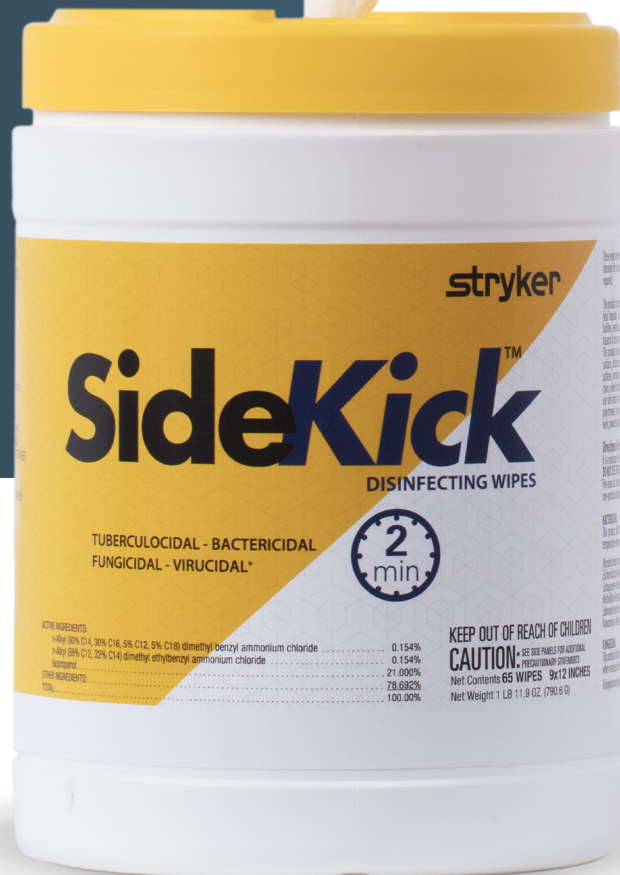

SideKick

Disinfecting Wipes

In the healthcare environment today, the high-touch surfaces surrounding our patients are at-risk areas for harboring viruses, bacteria, and other pathogens. Introducing SideKick, an easy-to-use wipe that is both convenient and highly effective at disinfecting pre-cleaned, hard, nonporous surfaces with its two-minute wet contact time. Its broad-spectrum efficacy is tuberculocidal, bactericidal, virucidal and fungicidal.



100%
compatible

We are your partners in prevention through the equipment we provide and in the solutions designed to take care of your products. SideKick Disinfecting Wipes are the only cleaning and disinfecting solution that is 100% compatible¹ for use on our beds, stretchers, patient support surfaces², furniture, chest compression systems, defibrillators, emergency care cots, cot fasteners and chairs.¹



2 minute
kill time

tuberculocidal, bactericidal,
virucidal and fungicidal

We've listened to your concerns and we understand your challenges. We have developed a product that meets your facility's many cleaning and disinfecting needs. **With SideKick wipes, you can feel confident standardizing with our products in hospital and EMS environments.**



Compatible

Tested 100% compatible on 200+ materials used in our beds, stretchers, patient support surfaces², furniture, chest compression systems, defibrillators, emergency care cots, cot fasteners and chairs.¹



Versatile

Broad-spectrum efficacy: tuberculocidal, bactericidal, fungicidal, virucidal.

Both an effective cleaner and disinfectant.²

Can be used on our products and other healthcare surfaces.

Neutral pH of 7.6.



Convenient

Two-minute wet contact time.

Ready-to-use premoistened wipe formulation: no mess or potential for chemical mixing errors.

SideKick Disinfecting Wipes

Specifications	Standard wipes	Large wipes
Model number	2060-000-001	2060-000-002
Wipe dimensions	6" x 10"	9" x 12"
Count per canister	100	65
Weight of canister	1 lb. 7.8 oz. (676.6g)	1 lb. 11.9 oz. (790.6g)

Standard features

- Compatible on Stryker's products¹
- Broad-spectrum efficacy
 - Tuberculocidal
 - Bactericidal
 - Fungicidal
 - Virucidal
- 2-minute wet contact time
- Neutral pH of 7.6
- Alcohol quat formulation
- For use on many healthcare surfaces²
- EPA reg number: 70144-2-92728

Active ingredients

- n-Alkyl (60% C14, 30% C16, 5% C12, 5% C18) dimethyl benzyl ammonium chloride: **0.154%**
- n-Alkyl (68% C12, 32% C14) dimethyl ethylbenzyl ammonium chloride: **0.154%**
- Isopropanol: **21.000%**
- Other ingredients: **78.692%**

Always refer to the SideKick wipes label for additional information and instructions for use. References and/or data for claims are on file at Stryker and available upon request.

1. Compatibility testing refers to the validation of Product surface materials performed by Stryker to not degrade, discolor or crack with exclusive use of SideKick for daily device cleaning and disinfection for the Product's stated expected life. Subject to footnote 2 below. Products refer to beds, stretchers, patient support surfaces, furniture, chest compression systems, defibrillators, emergency care cots, cot fasteners and chairs manufactured by and/or for Stryker.

2. Cleaning and disinfection claims exclude the Equilibrium support surface cover. Equilibrium by Dartex® and Dartex® are registered trademarks of Dartex Coatings, Ltd.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: SideKick, Stryker. All other trademarks are trademarks of their respective owners or holder.

3800 E. Centre Avenue
Portage, MI 49002 USA
t: 269 329 2100
toll free: 800 327 0770

stryker.com/sidekick