

# Sof.Care<sup>®</sup> Chair Cushion



## Convenient. Complete. Effective.

### 3-layer technology

Our Sof.Care family of products allows for a convenient, complete and effective way to help enhance pressure redistribution and patient comfort while helping improve hospital outcomes and prevent adverse events.

The Sof.Care Chair Cushion combines three layers of interlocking air cells to create a comfortably contoured support surface. Sof.Care elevates the body's sensitive pressure points off the surface of the chair and redistributes the patient's weight over the entire surface. The result is a superior combination of support and comfort.



### Sof.Care Chair Cushion

has air-channeling technology which redistributes a patient's weight to reduce tissue interface pressure.




### Sof.Care Bariatric Chair Cushion

has dual chambers for improved stability and optimal patient support. For nursing convenience, this cushion is also preinflated.

## Pressure management with 3-layer technology

The 3-layer air-channeling technology creates interlocking air cells to immerse and envelop the patient providing low tissue interface pressure. The 3-layer design provides total air support helping prevent the patient from bottoming out because the patient is always supported by a pocket of air.

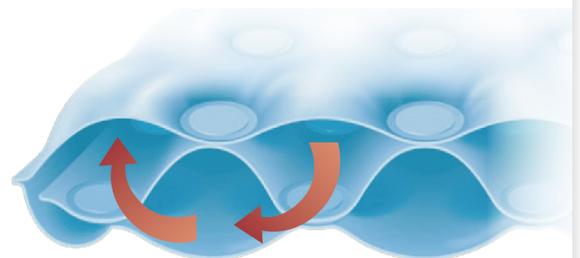

3-layer air-channeling technology

# Sof.Care Chair Cushion



## Sof.Care Chair Cushion features

- Air-channeling technology redistributes a patient's weight to reduce tissue interface pressure
- Protects sensitive bony prominences susceptible to pressure injuries
- Preinflated, ready to use
- Easy to transport
- Lightweight, waterproof and easy to clean
- Low-profile design fits all chairs
- Economical

## Sof.Care Bariatric Chair Cushion features

- Dual chambers for improved stability and optimal patient support
- Preinflated for nursing convenience
- Safe working load up to 700 lb.
- Air adjustment valves for customization if needed
- Small hand inflation pump available for pressure adjustments
- Same great features as the Sof.Care Chair Cushion
- Designed to fit most bariatric chairs

Specifications	Sof.Care Chair Cushion	Sof.Care Bariatric Chair Cushion
<b>Model number</b>	2770100000	2770200000
<b>Inflated length</b>	18" (47 cm)	24" (61 cm)
<b>Inflated width</b>	19.625" (50 cm)	32" (81 cm)
<b>Weight</b>	11.7 oz. (.33 kg)	2.25 lb. (1 kg)
<b>Weight capacity</b>	350 lb. (158.8 kg)	700 lb. (317.5 kg)
<b>Packaged</b>	6 per carton	5 per carton
<b>Warranty</b>	30 days	30 days
<b>Small hand inflation pump</b>	HCP700	HCP700

Stryker reserves the right to change specifications without notice. Single patient use only.

Stryker reserves the right to change specifications without notice. Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Sof.Care, Stryker. All other trademarks are trademarks of their respective owners or holder.

3800 E. Centre Avenue  
Portage, MI 49002 USA  
t: 269 329 2100  
toll free: 800 327 0770

[patientcare.stryker.com](http://patientcare.stryker.com)