

SPR Plus[®]

Low Air Loss Overlay




Convenient. Complete. Effective.

3-layer technology

Our Sof.Care family of products allows for a convenient way to assist with pressure redistribution. The SPR Plus overlay combines three layers of interlocking air cells to create a contoured support surface with low air loss microvents throughout.



Adjustable firmness setting allows caregiver to easily customize air support based on each patient's needs.

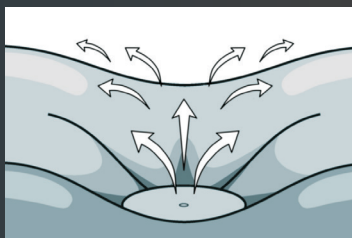


Designed with 12 direct low air loss microvents to help minimize moisture build up and help manage a patient's microclimate.



The low per-patient cost can reduce the need for expensive rental products.

Pressure management with 3-layer technology

The 3-layer air-channeling technology creates interlocking air cells to immerse and envelop the patient. The 3-layer design provides air support helping prevent the patient from bottoming out because the patient is supported by a pocket of air.

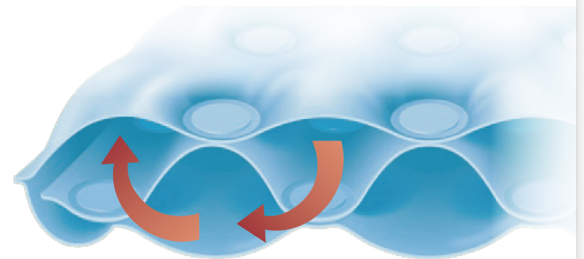

3-layer air-channeling technology

SPR Plus

Low Air Loss Overlay



SPR Plus Overlay features

- Air-channeling technology redistributes a patient's weight to reduce tissue interface pressure
- Helps protect sensitive bony prominences susceptible to pressure injuries
- Direct low air loss vents to help manage patient's microclimate of the skin
- Adjustable firmness settings
- Easy to transport
- Lightweight, waterproof and easy to clean
- Economical

Specifications

Specifications	SPR Plus Overlay
Model number	2790-100-000
Inflated length	78" (198 cm)
Inflated width	34" (86 cm)
Weight	6 lb. (2.72 kg)
Safe working load*	350 lb. (158.8 kg)
Flow rate	12.5 liters/minute
Warranty	30 days

Specifications

Specifications	Stryker Air Pump Plus
Model number	2863-000-000
Height	17.25" (43.8 cm)
Width	9.25" (23.5 cm)
Depth	7.25" (18.4 cm)
Weight	13.2 lb. (6.0 kg)

Single patient use only.

Ask about our other Sof.Care 3-layer air-channeling technology products. Stryker reserves the right to change specifications without notice.

Stryker reserves the right to change specifications without notice. Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Sof.Care, SPR Plus, Stryker, Stryker Air. All other trademarks are trademarks of their respective owners or holder.

3800 E. Centre Avenue
Portage, MI 49002 USA
t: 269 329 2100
toll free: 800 327 0770

stryker.com