

SPR Plus® Low Air Loss Overlay

3-layer Technology




Convenient. Complete. Effective.

Stryker's SofCare® family of products allows for a convenient, complete and effective way to assist pressure redistribution while improving hospital outcomes and helping to prevent facility acquired pressure ulcers.

The SPR Plus overlay combines three layers of interlocking air cells to create a contoured support surface with low air loss microvents throughout to help minimize moisture buildup and manage the patient microclimate.



Adjustable firmness setting allows caregiver to easily customize air support based on each patient's needs.

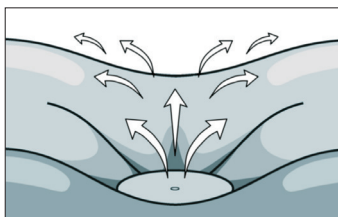


Designed with 12 direct low air loss microvents to help minimize moisture build up and manage a patient's microclimate.



The Low per-patient cost significantly reduces the need for expensive rental products.

Pressure Management with 3-layer Technology

The unique 3-layer air-channeling technology creates interlocking air cells to immerse and envelop the patient providing low tissue interface pressure. The 3-layer design provides total air support preventing the patient from bottoming out because the patient is always supported by a pocket of air.

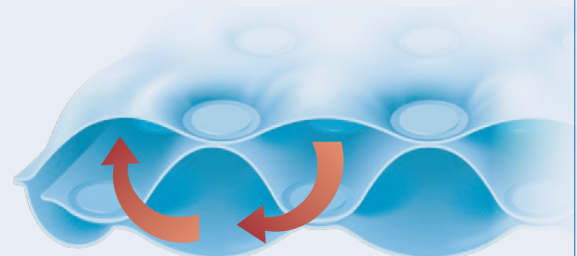

3-Layer Air-Channeling Technology

SPR Plus Overlay



SPR Plus Overlay Features

- Unique air-channelling technology redistributes a patient's weight to reduce tissue interface pressure
- Protects sensitive bony prominences susceptible to pressure ulcers
- Direct low air loss vents to help manage patient's microclimate
- Adjustable firmness settings
- Easy to transport
- Lightweight, waterproof and easy to clean
- Economical

Specifications

SPR Plus Overlay

Model Number	2790-100-100
Inflated Length	78" (198 cm)
Inflated Width	34" (86 cm)
Weight	6 lb. (2.72 kg)
Weight Capacity	350 lb. (158.8 kg)
Flow Rate	12.5 liters/minute
Warranty	30 days

Control Unit

Stryker Air[™] Pump

Model Number	2861-000-026
Width	10" (25.4 cm)
Length	10" (25.4 cm)
Weight	5.5 lb (2.5 kg)

Single Patient Use Only

Ask about our other SofCare 3-layer air-channelling technology products. Stryker reserves the right to change specifications without notice.

3800 E. Centre Ave.
Portage, MI 49002 U.S.A.
t: 269 389 2100 f: 269 329 2311
toll free: 800 787 9537

www.stryker.com