

**stryker**

# Maximizing minimally invasive surgery

Signature Portfolio  
MIS Attachments



# Maximizing MIS performance

Your needs change from day to day. That's precisely why our Signature Portfolio high speed drill platform was designed to offer you the depth and breadth you require in minimally invasive surgery (MIS) attachments and cutting accessories. Developed from user input, our MIS attachments have been rebalanced and refined to provide you with a hands-on experience that's completely customizable to your preferences.

To maximize your performance even further, we've optimized and expanded our offering of MIS cutting accessories to give you access to the largest selection of MIS cutting accessories – including the largest MIS cutting accessory head on the market.

Each MIS cutting accessory has been designed to provide multiple options for the precise reach, customization and exposure you need for a particular patient or procedure. Large diameter cutting accessories offer four different telescoping exposure points (in 1mm increments) and small diameter cutting accessories offer ten; all easily adjusted using built-in notches for exposure and minimum/maximum lengths.

This flexibility to extend and retract each cutting tool allows you to cost effectively secure the broad options and functionality you want through a streamlined set of instruments in your inventory or surgical field.


**The only MIS cutting accessories on the market with multi-notch (4-10mm) telescoping ability; maximizes length options while minimizing cost and inventory**



◀ Patented 2Flute cutting accessory, smooth and precise for the rapid dissection of bone



◀ Head diameter up to 5mm, largest on the market

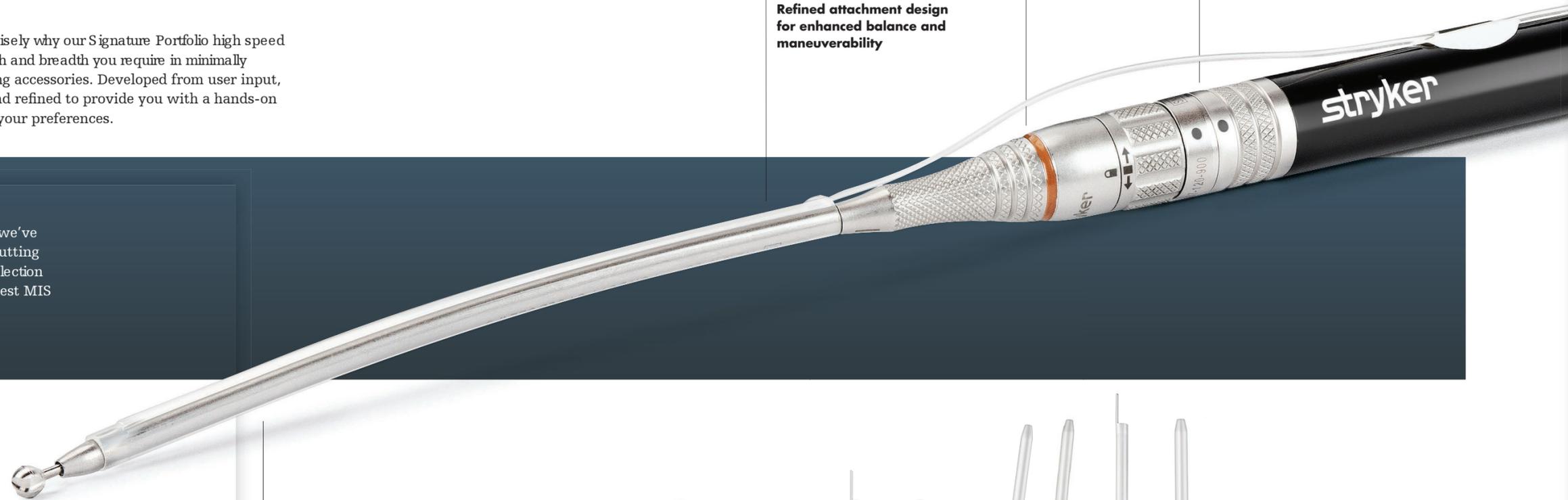

Longer, concave design for enhanced line of sight and fine-edge knurling for better grip and more options for positioning your hand

Refined attachment design for enhanced balance and maneuverability

Color coding instantly identifies matching attachments and cutting accessories

Alignment dots and Snap-Connect guide assembly/disassembly of attachments to motor in an easy and consistent manner

## NEW! Irrigation Sleeves

Deliver precise irrigation by simply sliding over your 13cm or 16cm MIS attachment for a total MIS solution

Low profile with ergonomic design facilitates site access and visibility

Sterile, disposable and compatible with CORE Console disposable irrigation cassettes



## Portfolio options

Signature Portfolio offers more than 50 MIS cutting accessories — including two recent additions to our lineup — in an unparalleled variety of head types, diameters, (1.5 - 5.0mm) and lengths (13-16cm). For a complete listing, see our back cover or contact your sales representative.

Match Head



Precision Match Head



Diamond Match Head



Coarse Diamond Match Head



Round



Precision Round



Coarse Diamond




Diamond Extra Coarse





## The perfect touch

When used with  $\pi$ drive or  $\pi$ drive+ motors, your cutting accessories can be optimized even further using I.D. Touch Software on the CORE 2 Console instrument driver. This industry exclusive software enables you to adjust torque from 0-100% to customize how your electric drill responds when pressure is applied during use — from low to high sensitivity — redefining performance to your liking. Your NSE sales representative can guide you through the options and save the performance and feel that's right for you in user preferences.

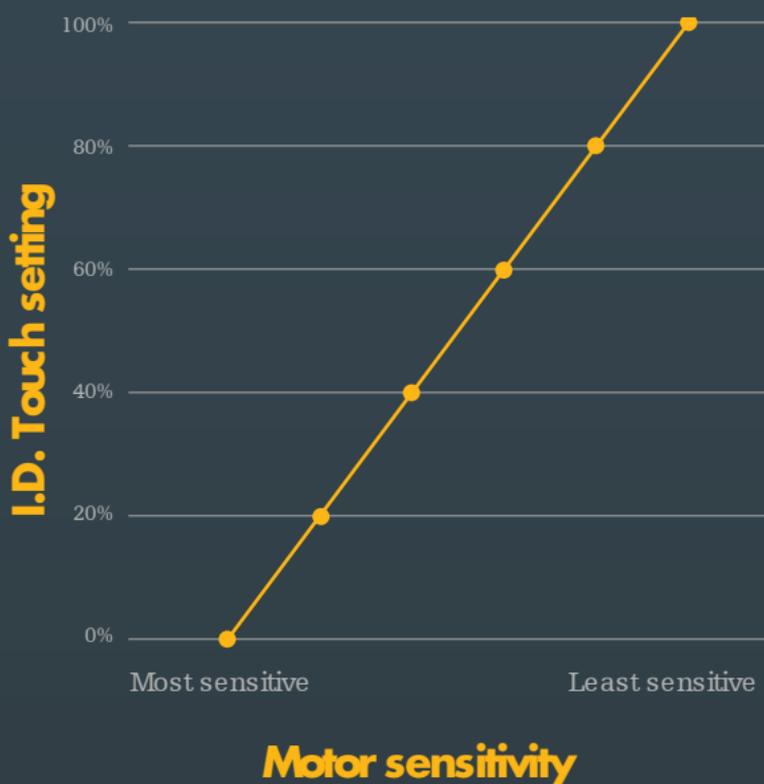

## Signature Portfolio MIS Attachments and Cutting Accessories

To order our MIS attachments, cutting accessories or Irrigation Sleeves, call 800 253 3210.

### Attachments

Part number	Product description	Color band
5407-120-900	MIS Hub	
5407-120-900A	MIS Angled Hub	
5407-120-950	13cm MIS Straight Nose Tube	Lilac
5407-120-950A	13cm MIS Angled Nose Tube	Lilac
5407-120-950C	13cm MIS Curved Nose Tube	Lilac
5407-120-950F	13cm MIS Modified Curve Nose Tube	Lilac
5407-120-950H	13cm MIS Straight Hooded Nose Tube	Lilac
5407-120-970	16cm MIS Straight Nose Tube	Brown
5407-120-970A	16cm MIS Angled Nose Tube	Brown
5407-120-970C	16cm MIS Curved Nose Tube	Brown
5407-120-970H	16cm MIS Straight Hooded Nose Tube	Brown

### Cutting Accessories

8450 (lilac) series = 13cm shank      8470 (brown) series = 16cm shank

Part number	Product description	Color band
8450-107-015	1.5mm Match Head	Lilac
8470-107-015	1.5mm Match Head	Brown
8450-107-020	2.0mm Match Head	Lilac
8470-107-020	2.0mm Match Head	Brown
8450-107-025	2.5mm Match Head	Lilac
8470-107-025	2.5mm Match Head	Brown
8450-107-025L	2.5mm Extended Match Head	Lilac
8450-107-030	3.0mm Match Head	Lilac
8470-107-030	3.0mm Match Head	Brown
8450-107-525	2.5mm Precision Match Head	Lilac
8470-107-525	2.5mm Precision Match Head	Brown
8450-107-525L	2.5mm Extended Precision Match Head	Lilac
8450-107-530	3.0mm Precision Match Head	Lilac
8470-107-530	3.0mm Precision Match Head	Brown
8450-107-015D	1.5mm Diamond Match Head	Lilac
8470-107-015D	1.5mm Diamond Match Head	Brown
8450-107-020D	2.0mm Diamond Match Head	Lilac
8470-107-020D	2.0mm Diamond Match Head	Brown
8450-107-025D	2.5mm Diamond Match Head	Lilac
8470-107-025D	2.5mm Diamond Match Head	Brown
8450-107-030D	3.0mm Diamond Match Head	Lilac
8470-107-030D	3.0mm Diamond Match Head	Brown

### Cutting Accessories continued

Part number	Product description	Color band
8450-107-035D	3.5mm Diamond Match Head	Lilac
8470-107-035D	3.5mm Diamond Match Head	Brown
8450-107-040D	4.0mm Diamond Match Head	Lilac
8470-107-040D	4.0mm Diamond Match Head	Brown
8450-107-030DC	3.0mm Coarse Diamond Match Head	Lilac
8470-107-030DC	3.0mm Coarse Diamond Match Head	Brown
8450-107-040DC	4.0mm Coarse Diamond Match Head	Lilac
8470-107-040DC	4.0mm Coarse Diamond Match Head	Brown
8450-010-020	2.0mm Round	Lilac
8470-010-020	2.0mm Round	Brown
8450-010-030	3.0mm Round	Lilac
8470-010-030	3.0mm Round	Brown
8450-010-040	4.0mm Round	Lilac
8470-010-040	4.0mm Round	Brown
8450-009-020	2.0mm Precision Round	Lilac
8470-009-020	2.0mm Precision Round	Brown
8450-009-030	3.0mm Precision Round	Lilac
8470-009-030	3.0mm Precision Round	Brown
8450-009-040	4.0mm Precision Round	Lilac
8470-009-040	4.0mm Precision Round	Brown
8450-009-050	5.0mm Precision Round	Lilac
8470-009-050	5.0mm Precision Round	Brown
8450-012-020DC	2.0mm Coarse Diamond	Lilac
8470-012-020DC	2.0mm Coarse Diamond	Brown
8450-012-030DC	3.0mm Coarse Diamond	Lilac
8470-012-030DC	3.0mm Coarse Diamond	Brown
8450-012-040DC	4.0mm Coarse Diamond	Lilac
8470-012-040DC	4.0mm Coarse Diamond	Brown
8450-013-040DX	4.0mm Diamond Extra Coarse	Lilac
8470-013-040DX	4.0mm Diamond Extra Coarse	Brown
8450-013-050DX	5.0mm Diamond Extra Coarse	Lilac
8470-013-050DX	5.0mm Diamond Extra Coarse	Brown

### Irrigation Sleeves

Part number	Product description
5407-010-950	For use with 13cm (lilac) MIS nose tubes
5407-010-970	For use with 16cm (brown) MIS nose tubes

## Neuro Spine ENT

This document is intended solely for the use of healthcare professionals. A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. We do not dispense medical advice and recommend that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate Stryker's products. A surgeon must always refer to the package insert, product label and/or instructions for use, including the instructions for cleaning and sterilization (if applicable), before using any of Stryker's products. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your representative if you have questions about the availability of Stryker's products in your area.

Stryker Corporation or its affiliates own, use, or have applied for the following trademarks or service marks: CORE, I.D. Touch,  $\pi$ drive and Stryker. All other trademarks are trademarks of their respective owners or holders.

Stryker  
4100 East Milham Avenue  
Kalamazoo, MI 49001 USA  
t: 269 323 7700 f: 269 323 2742  
toll free: 800 253 3210  
nse.stryker.com

D0000002511AB  
G55/PS  
Copyright © 2018 Stryker  
Printed in USA