

stryker



Driving
performance,

**delivering
versatility**

**Spetzler-Malis and SILVERGlide
Non-Stick Bipolar Forceps**

Your **superior experience** starts with our people

Supporting you throughout your procedures

People

Our people stand out. They're dedicated, highly skilled experts that care deeply and listen carefully. When you need a partner in improving care, we're by your side.

Culture

We create a dynamic, supportive environment that enables amazing people to do amazing things. Count on us to do what's right and deliver—every time.

Product innovation

Innovation is in our DNA—from Dr. Homer Stryker to the researchers, engineers and surgeons creating solutions that are reshaping neurotechnology today.

Medical education

Innovation and skill come together in our educational outreach to providers. Each year, we help surgical teams around the world master the latest medical technologies and care practices.

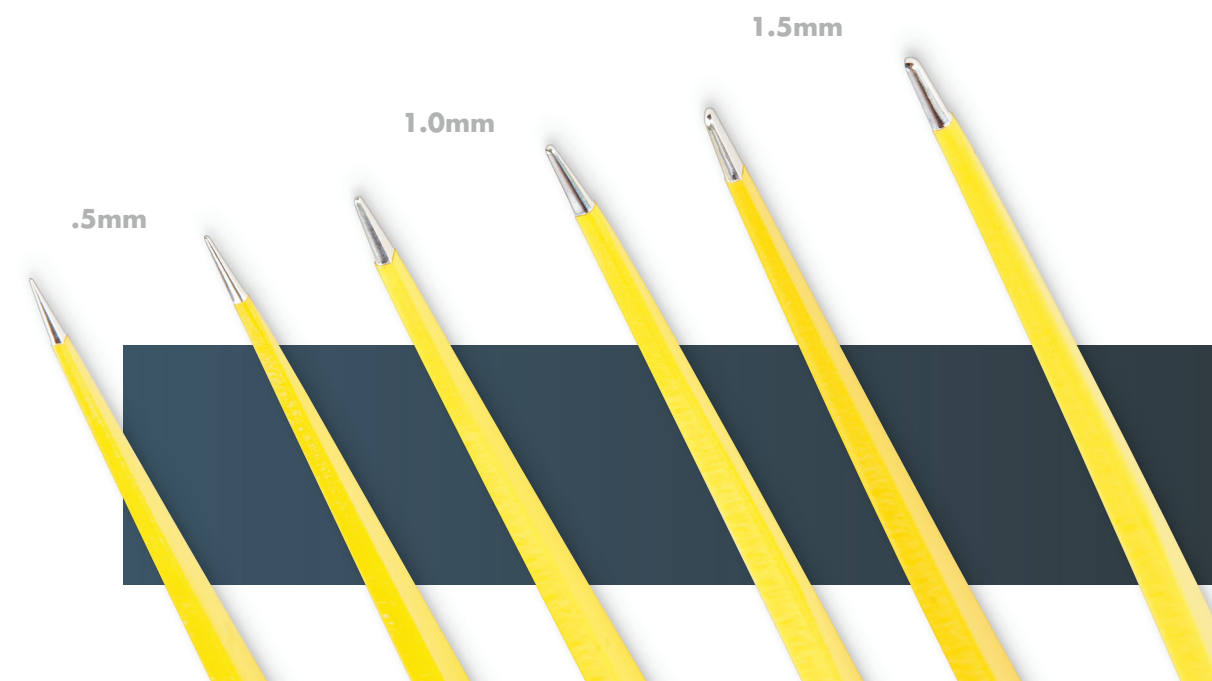
Research and development

We amplify our investments in R&D by listening to our customers and collaborating with them to create solutions that are simply superior.

Outcomes

We share your values and your commitment to innovative care that heals effectively and efficiently. With you, we're driven to make healthcare better.

All of our non-stick bipolar forceps **resist tissue buildup** and provide **excellent visualisation**. Regardless of your preference for single-use Spetzler-Malis or reusable SILVERGlide, you can expect **precision, control** and a responsive extension of your hands and eyes from us.



Disposable Spetzler-Malis

Standard Bipolar Forceps

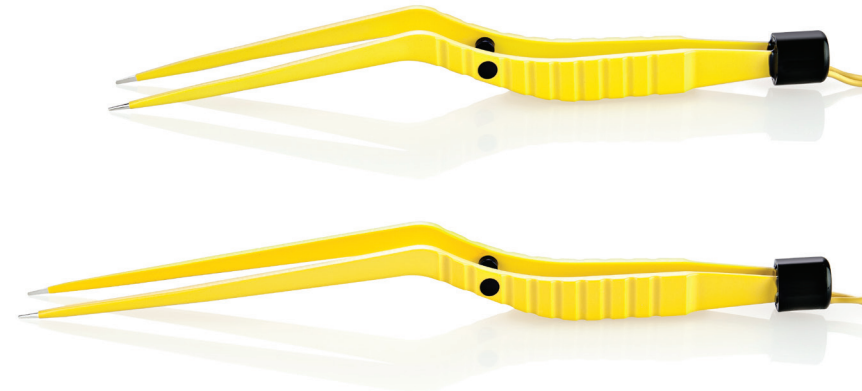
Flexible lightweight cord and remote connection plug helps achieve balance

Spring force elevates control and tactile feedback

Single use ensures excellent tip alignment and coagulation performance

Sleek bayonet design improves visualisation

Non-stick technology enhances performance



7" (18cm)

6760-180-005 0.5mm tip
6760-180-010 1.0mm tip
6760-180-015 1.5mm tip

8" (20cm)

6760-200-005 0.5mm tip
6760-200-010 1.0mm tip
6760-200-015 1.5mm tip

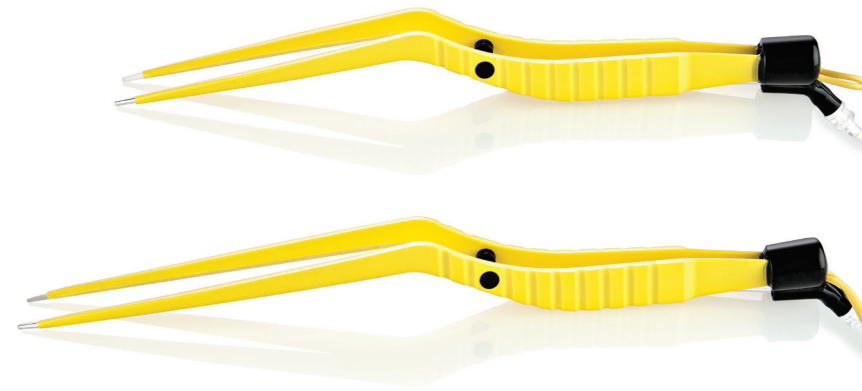
9" (23cm)

6760-230-005 0.5mm tip
6760-230-010 1.0mm tip
6760-230-015 1.5mm tip

Dual Irrigating Bipolar Forceps

Transparent fine tubing allows for better observation

Dual irrigation design maintains moisture of the tips and vessels



7" (18cm)

6770-180-005 0.5mm tip
6770-180-010 1.0mm tip
6770-180-015 1.5mm tip

8" (20cm)

6770-200-005 0.5mm tip
6770-200-010 1.0mm tip
6770-200-015 1.5mm tip

9" (23cm)

6770-230-005 0.5mm tip
6770-230-010 1.0mm tip
6770-230-015 1.5mm tip

Slim Bipolar Forceps

Slim profile provides access in narrow fields



8" (20cm)

6780-200-005 0.5mm tip
6780-200-010 1.0mm tip
6780-200-015 1.5mm tip

9" (23cm)

6780-230-005 0.5mm tip
6780-230-010 1.0mm tip
6780-230-015 1.5mm tip

SILVERGlide Non-Stick

ProSeries Bipolar Forceps



Standard two-pin round connector

works with most standard foot-switching bipolar electro-surgical generators

Silver plating
reduces collateral heating, "hot spots" and sticking for performance

Ultra-thin tines
allow better access in surgery

Lightweight handle
for easy manipulation and balance

Guide-stop
minimises scissoring and increases control



5.7" (14.5cm) bayonet

- 6730-145-004 0.4mm
- 6730-145-007 0.7mm
- 6730-145-010 1.0mm
- 6730-145-012 1.2mm
- 6730-145-015 1.5mm

6.3" (16cm) bayonet

- 6730-160-004 0.4mm
- 6730-160-007 0.7mm
- 6730-160-010 1.0mm
- 6730-160-012 1.2mm
- 6730-160-504, 45° up 0.4mm

7.5" (19cm) bayonet

- 6730-190-004 0.4mm
- 6730-190-007 0.7mm
- 6730-190-010 1.0mm
- 6730-190-012 1.2mm
- 6730-190-015 1.5mm
- 6730-190-020 2.0mm
- 6730-190-504, 45° up 0.4mm
- 6730-190-512, 45° up 1.2mm
- 6730-190-612, 45° down 1.2mm

8.6" (22cm) bayonet

- 6730-220-004 0.4mm
- 6730-220-007 0.7mm
- 6730-220-012 1.2mm
- 6730-220-015 1.5mm
- 6730-220-020 2.0mm
- 6730-220-107, 15° up 0.7mm
- 6730-220-307, 30° up 0.7mm
- 6730-220-507, 45° up 0.7mm
- 6730-220-607, 45° down 0.7mm
- 6730-220-512, 45° up 1.2mm
- 6730-220-612, 45° down 1.2mm

7.5" (19cm) straight bayonet

- 6735-190-007 0.7mm
- 6735-190-010 1.0mm
- 6735-190-310, 30° up 1.0mm
- 6735-190-320, 30° up 2.0mm



9.4" (24cm) bayonet

- 6730-240-004 0.4mm
- 6730-240-007 0.7mm
- 6730-240-012 1.2mm
- 6730-240-015 1.5mm
- 6730-240-020 2.0mm
- 6730-240-507, 45° up 0.7mm
- 6730-240-607, 45° down 0.7mm
- 6730-240-807, 80° up 0.7mm
- 6730-240-907, 80° down 0.7mm

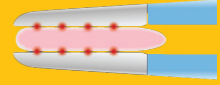
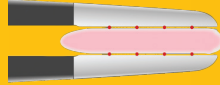
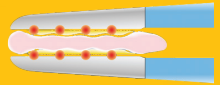
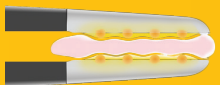
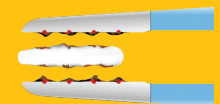



11" (28cm) bayonet

- 6730-280-012 1.2mm
- 6730-280-015 1.5mm

What sets SILVERGlide apart?

The elimination of tissue buildup allows you to work without interruption. These non-stick characteristics continue to perform use after use. The SILVERGlide alloy's level of thermal diffusivity is 37 times higher than that of stainless steel commonly used in forceps.

| | Standard forceps | SILVERGlide Non-Stick | |
|--|---|---|--|
| "Hot spots" from sparks occur during coagulation |  |  | "Hot spots" start to occur during coagulation |
| The "hot spots" cause localised sticking |  |  | Heat is diffused rapidly to prevent sticking |
| Tissue seared to the forceps tips is ripped away |  |  | The forceps can be removed without tissue sticking |

Accessories

SILVERGlide Bipolar Forceps

SILVERGlide Bipolar Forceps sterilisation trays

6700-002-000 Two forceps sterilisation tray
6700-005-000 Five forceps sterilisation tray
6700-010-000 10 forceps sterilisation tray (forceps not included)

Other SILVERGlide accessories

6700-000-002 Polishing cloths
6700-000-001 Tip guards
6701-000-000 Reusable cable



Disposable Spetzler-Malis Bipolar Forceps cord and irrigation sets

Rotary cord and tubing sets

6790-100-003 Bipolar irrigator integrated tubing and cord set, disposable
6790-100-004 Bipolar irrigator integrated tubing and cord set, disposable, unitised plug

Gravity style cord and tubing sets

6790-100-001 Integrated tubing and bipolar cord set, disposable
6790-100-002 Integrated tubing and bipolar cord set, disposable, unitised plug

Cord only sets

6790-100-005 Bipolar cord set, disposable
6790-100-006 Bipolar cord set, disposable, unitised plug



Malis Bipolar Generator and Irrigator system

Malis Bipolar Generator and Irrigator system

6700-313-500 Malis Bipolar Generator (includes power cord, does not include foot pedal)
6700-313-600 Malis Bipolar Irrigator (does not include power cord)
0996-851-011 Universal power cord
6700-313-004 Bipolar generator/irrigator interconnect cable
6700-313-700 Bipolar generator cart



Irrigator
Optional integrated irrigator

Oversized color screen
for easy comprehension

Programmable user profiles
store information for up to eight surgeons

CUT and COAG knobs
facilitate quick setting adjustments

Independent bipolar outputs
enhance functionality

DUAL WAVE technology
provides confidence during cutting and coagulating

Foot pedals

6700-313-003 Single foot pedal (COAG only, does not include required adapter cable)
6700-313-304 Single foot pedal adapter cable
6700-313-002 Dual foot pedal



**Whether you choose reusable or disposable,
our non-stick bipolar forceps act as a specialised
extension of your hands and eyes.**

For more information, contact your Neurosurgical Sales Representative.

Neurosurgical

This document is intended solely for the use of healthcare professionals. A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. We do not dispense medical advice and recommend that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate Stryker's products. A surgeon must always refer to the package insert, product label and/or instructions for use, including the instructions for cleaning and sterilisation (if applicable), before using any of Stryker's products. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your representative if you have questions about the availability of Stryker's products in your area.

Stryker or its affiliated entities own, use, or have applied for the following trademarks or service marks: Malis, SILVERGlide, Spetzler and Stryker. All other trademarks are trademarks of their respective owners or holders.

D0000053188 Rev. AA_SSP_AUS
G55/PS
Copyright © 2020 Stryker
Printed in USA

Stryker Australia
8 Herbert Street
St.Leonards
NSW 2065 Australia
Ph: +61 2 9467 1000
Fax: +61 2 9467 1010

Stryker New Zealand
515 Mt. Wellington Highway
Auckland 1060 New Zealand
Ph: +64 9 573 1890
Fax: +64 9 573 1891