


Everest[®]

Degenerative Spinal System

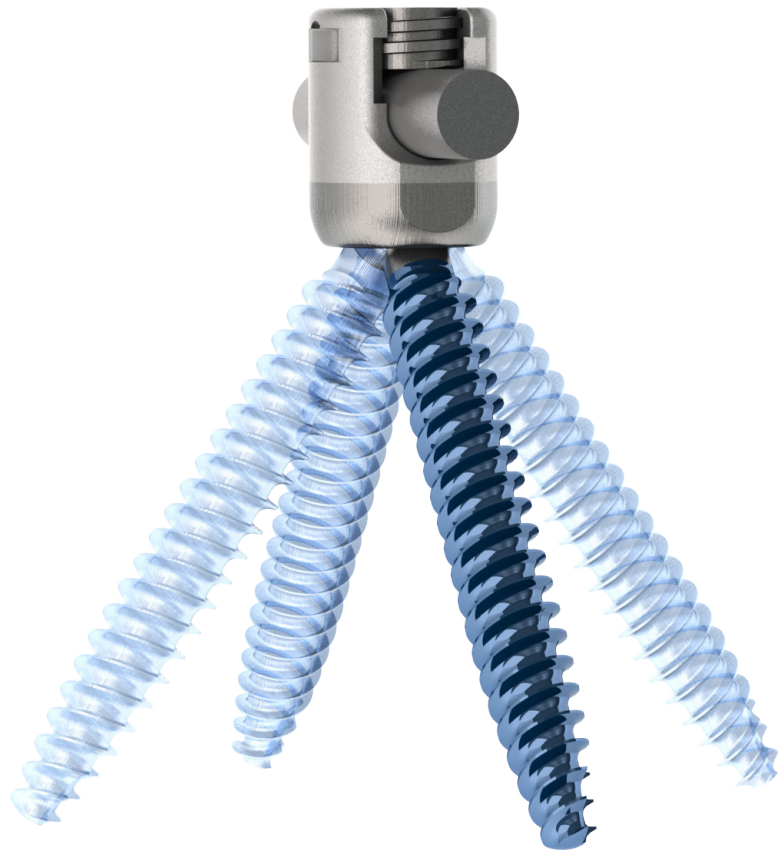

The Everest Degenerative Spinal System is a top-loading polyaxial pedicle screw system with the ability to accommodate titanium and cobalt chrome rods of two different diameters. The dual-lead thread pattern allows for faster insertion and demonstrates increased pullout strength.¹ The mixed-metal (Ti/CoCr) tulip minimized head splay and demonstrated improved mechanical performance when tested against an all-titanium alloy screw.^{2,3,4}

1. TR-486. Data on file.
2. TR-595. Data on file.
3. TR-566. Data on file.
4. TR-576. Data on file.

Everest Degenerative Spinal System

Variable rod compatibility



- Screw head accepts both Ø5.5 and Ø6.0mm rods
- Set contains pre-lordosed rods in 5mm increments

Mixed material



- Mixed-metal (Ti/CoCr) tulip minimized head splay and demonstrated improved mechanical performance when tested against an all-titanium alloy screw^{2,3,4}
- Provides a posable friction head

Simplified Set Screw



- Everest Set Screw features a modified square thread design, facilitating set screw introduction
- 9° allowable misalignment angle

Optimized thread pattern



- Dual-lead thread pattern for faster insertion and increased pullout strength¹

1. TR-486. Data on file.
2. TR-595. Data on file.
3. TR-566. Data on file.
4. TR-576. Data on file.

Spine division

This document is intended solely for the use of healthcare professionals. A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A surgeon must always refer to the package insert, product label and/or instructions for use before using any Stryker product. The products depicted are CE marked according to the Medical Device Directive 93/42/EEC. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Everest, K2M, Stryker. All other trademarks are trademarks of their respective owners or holders.

