


# Mesa®

## Small Stature Spinal System

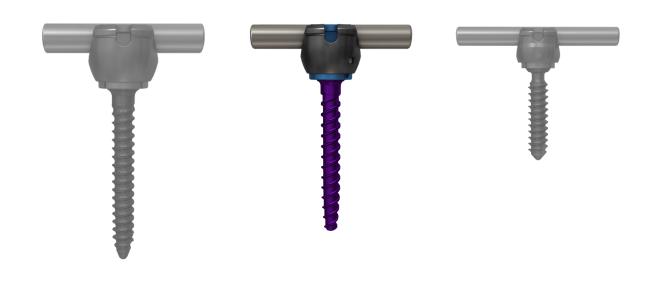

The Mesa Small Stature Spinal System includes a top-loading, low-profile, next generation Mesa screw featuring Zero-Torque Technology, coupled with light ergonomic instrumentation. This system is poised to address difficult correction maneuvers in small stature patients.

### Mesa Small Stature Spinal System

#### **Implant details**



The Mesa Small Stature Screw is 10% slimmer than the Mesa Deformity Screw



The Mesa Small Stature Screw is 24% lower than the Mesa Deformity Screw



The Mesa Small Stature Screw Head is 18% smaller than the Mesa Deformity Screw Head

#### **Reduction options**



Deformity Cricket The Deformity Reduction Jack (Cricket) is designed to effectively eliminate the need for reduction screws by providing 30mm of travel distance to the rod and allowing for simultaneous multi-axial translation and reduction.



The single action rod reducer (Firefly) offers 10mm of reduction in a pistol grip instrument design.



Dragonfly

The Dual Action Rod Reducer (Dragonfly) offers 15mm of single action reduction and also partially locks the Mesa screw.

#### Ability to segmentally or globally derotate spine to achieve axial plane correction

- Deformity Reduction Jack (Cricket) provides ability to accomplish correction maneuvers in all planes
- 65° range of motion
- Ergonomic instrumentation is designed to allow for slow, controlled correction of spine, while distributing forces across entire construct

#### **Transition options**

Complete offering of transition options allowing connection to  $\emptyset 3.5$ , 4.5, and 5.5mm systems.







Parallel Connectors



Transition Rods

#### Spine division

This document is intended solely for the use of healthcare professionals. A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A surgeon must always refer to the package insert, product label and/or instructions for use before using any Stryker product. The products depicted are CE marked according to the Medical Device Directive 93/42/EEC. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Cricket, K2M, Mesa, Stryker, Zero-torque Technology. All other trademarks are trademarks of their respective owners or holders.



K2M, Inc. 600 Hope Parkway SE Leesburg, VA 20175 t: 571 919 2000

www.stryker.com