


stryker

AlphaVent Knotless

Self-punching, knotless anchor platform

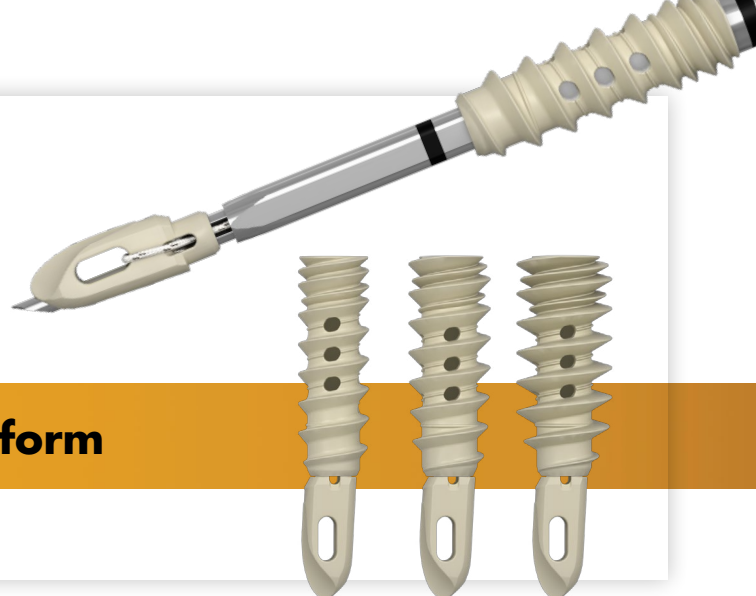

4.75mm 5.5mm 6.5mm

Streamlined design




Coupled, self-punching anchor system

Alleviates need for pilot hole creation*

Coupled eyelet and anchor facilitate streamlined tensioning


*Pilot hole creation required for hard bone indications.

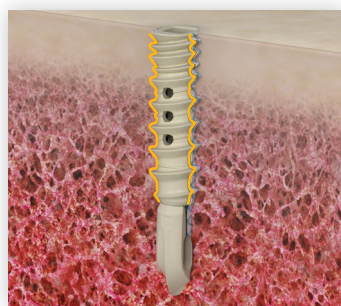

Dual-thread Technology

Designed to support anchor performance through cortical thread engagement¹

Provides tactile feedback during cortical thread engagement¹

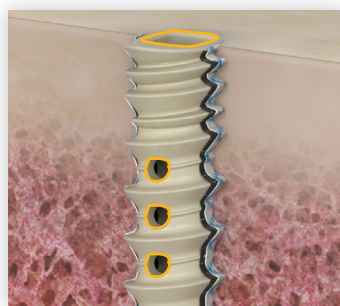

Complete venting and cannulation

Allows exchange of blood and other autologous cells at the implantation site




Your repair. Your way.

Accommodates a variety of suture configurations

Anchor length enables use in a single and double row repair

Features

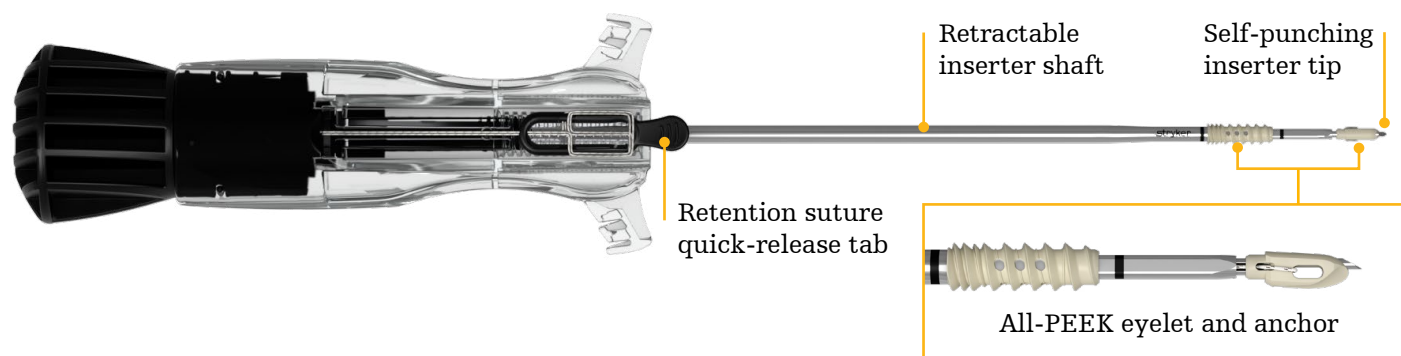

Hard bone instrumentation



3.9mm Omega and 4.75mm AlphaVent awl



4.75mm AlphaVent and Omega awl



4.75mm AlphaVent punch-tap



4.75mm AlphaVent tap



5.5mm AlphaVent tap



4.75mm AlphaVent and Omega drill



5.5mm AlphaVent and 4.75mm hard bone Omega drill

AlphaVent Knotless offerings

Implants

3911956547	AlphaVent Knotless SP, 4.75mm PEEK anchor
3911956555	AlphaVent Knotless SP, 5.5mm PEEK anchor
3911956565	AlphaVent Knotless SP, 6.5mm PEEK anchor

Instrumentation

3910947200	3.9mm Omega and 4.75mm AlphaVent awl
3910947201	4.75mm AlphaVent and Omega awl
3910947202	4.75mm AlphaVent and Omega drill
3910947207	4.75mm AlphaVent tap
3910947208	4.75mm AlphaVent punch-tap
3910947203	5.5mm AlphaVent and 4.75mm hard bone Omega drill
3910947209	5.5mm AlphaVent tap

References:

1. Data on file TR22963

Sports Medicine

This document is intended solely for the use of healthcare professionals. A healthcare professional must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that healthcare professionals be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A healthcare professional must always refer to the package insert, product label and/or instructions for use before using any Stryker product. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: AlphaVent, Omega and Stryker. All other trademarks are trademarks of their respective owners or holders.

The products depicted are CE marked in accordance with applicable EU Regulations and Directives.

SM-SHI-BROC-2093074

Copyright © 2025 Stryker
stryker.com

Stryker European Operations Ltd.
Herikerbergweg 110
1101 CM Amsterdam
The Netherlands

Manufactured by:
Stryker Endoscopy
5900 Optical Court
San Jose, CA 95138 USA
1-800-624-4422

CE Mark/Notified Body:
CE0197
TUV Rheinland