

# XBraid<sub>TT</sub>

## Overview

- XBraid<sub>TT</sub> is a non-absorbable, smooth, low profile, 100% ultra-high molecular weight polyethylene suture tape.
- Designed to increase tendon-to-bone contact especially where poor tissue quality may be a concern.
- XBraid<sub>TT</sub> has a wide suture footprint allowing for uniform pressure along the surface area of the repair.
- Tie-able tape with comparable knot stacks to standard #2 suture. XBraid<sub>TT</sub> easily slides after anchor insertion and during knot tying.
- Easy to open packaging designed for efficiency in the operating room.
- Custom braid technology provides a suture tape that is specifically designed to not fold under tension.
- XBraid<sub>TT</sub> will not pull through tissue as easily as #2 round suture<sup>1</sup>.
  - XBraid<sub>TT</sub> ultimate tissue pull through load: 206N<sup>1</sup>
  - #2 suture ultimate tissue pull through load: 141N<sup>1</sup>




## Dimensions

- 2.0mm XBraid<sub>TT</sub> Blue Co-Braid: 39" Strand
- 1.2mm XBraid<sub>TT</sub> Blue Co-Braid: 45" Strand

## Tips and Instructions

- Packaging
  1. Open packaging from chevron end.
  2. Remove the suture carrier card.
  3. Pull the exposed suture loop from the packaging card in order to get the suture out.

Part Number	Description
3910-900-018	2.0mm XBraidTT, Blue Co-Braid 39" Strand
3910-900-017	1.2mm XBraidTT, Blue Co-Braid 45" Strand

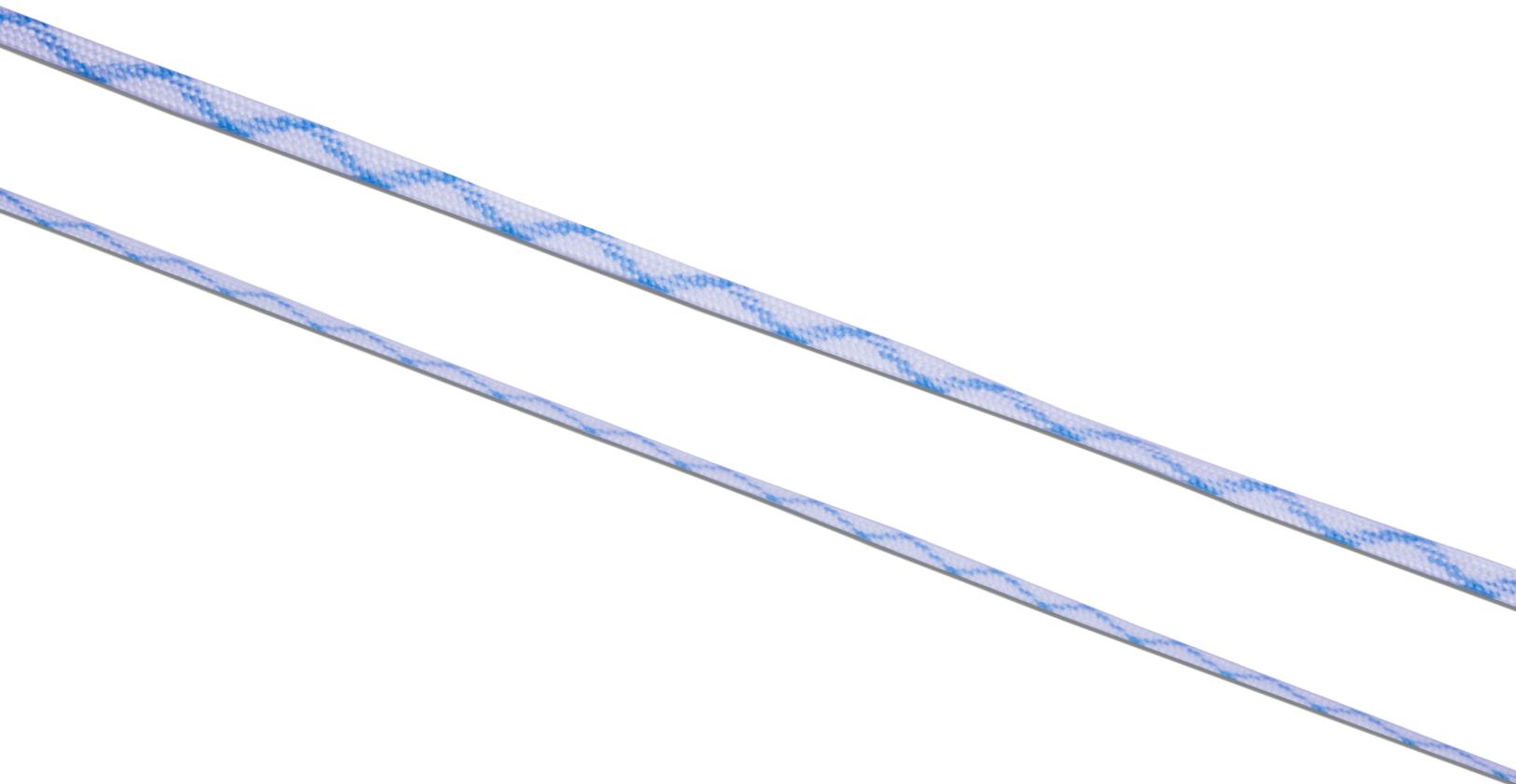

References:

1. Test Report: DHFD12758

---

5670 Greenwood Plaza Blvd. Ste. 200  
Greenwood Village, CO 80111  
t: 201 831 5000  
www.stryker.com  
www.sportsmedicine.stryker.com

A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A surgeon must always refer to the package insert, product label and/or instructions for use before using any Stryker product. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks; XBraid TT, Stryker. All other trademarks are trademarks of their respective owners or holders.

1000901974 Rev. A  
Copyright © 2015 Stryker  
Printed in USA